

An optimized comet-based in vitro DNA repair assay to assess base and nucleotide excision repair activity

Peer-reviewed author version

Vodenkova, Sona; Azqueta, Amaya; Collins, Andrew; Dusinska, Maria; Gaivao, Isabel; Moller, Peter; Opattova, Alena; Vodicka, Pavel; Godschalk, Roger W. L. & LANGIE, Sabine (2020) An optimized comet-based in vitro DNA repair assay to assess base and nucleotide excision repair activity. In: Nature protocols, 15 (12) , p. 3844 -3878.

DOI: 10.1038/s41596-020-0401-x

Handle: <http://hdl.handle.net/1942/33016>

TITLE: The comet-based *in vitro* DNA repair assay: a standardized method to assess an individual's DNA repair activity.

Authors:

Sona Vodenkova*, Department of Molecular Biology of Cancer, Institute of Experimental Medicine of the Czech Academy of Sciences, Prague; and Department of Medical Genetics, Third Faculty of Medicine, Charles University, Prague, Czech Republic, sona.vodenkova@iem.cas.cz,

Amaya Azqueta Osoz*, Department of Pharmacology and Toxicology, University of Navarra, and IdiSNA, Navarra Institute for Health Research, Pamplona, Spain, amazqueta@unav.es,

Andrew Collins, Department of Nutrition, University of Oslo, Norway, a.r.collins@medisin.uio.no,

Maria Dusinska, Department of Environmental Chemistry, Health Effects Laboratory, NILU-Norwegian Institute for Air Research, Kjeller, Norway, mdu@nilu.no,

Isabel O'Neill De Mascarenhas Gaivão, Genetics and Biotechnology Department and Veterinary and Animal Research Centre (CECAV), Universidade de Trás-os-Montes e Alto Douro, Vila Real, Portugal, igaivao@utad.pt,

Peter Møller, Department of Public Health, Section of Environmental Health, University of Copenhagen, Copenhagen, Denmark, pemo@sund.ku.dk,

Alena Opattova, Department of Molecular Biology of Cancer, Institute of Experimental Medicine of the Czech Academy of Sciences, Prague, Czech Republic, alena.opattova@iem.cas.cz,

Pavel Vodicka, Department of Molecular Biology of Cancer, Institute of Experimental Medicine of the Czech Academy of Sciences, Prague; and Biomedical Center, Medical faculty in Pilsen, Charles University in Prague, Czech Republic, pavel.vodicka@iem.cas.cz,

Roger W.L. Godschalk, Department of Pharmacology & Toxicology, School for Nutrition and Translational Research in Metabolism (NUTRIM), Maastricht University, The Netherlands, r.godschalk@maastrichtuniversity.nl,

Sabine A.S. Langie[#], Centre for Environmental Sciences, Hasselt University, Hasselt, Belgium, sabine.langie@uhasselt.be; and Department of Pharmacology & Toxicology, School for Nutrition and Translational Research in Metabolism (NUTRIM), Maastricht University, The Netherlands, s.langie@maastrichtuniversity.nl.

* shared first author

[#]Corresponding author

Please address all correspondence to:

Sabine Langie

Department of Pharmacology & Toxicology

School for Nutrition and Translational Research in Metabolism (NUTRIM)

Maastricht University, Universiteitssingel 50, PO Box 616

6200MD, Maastricht, The Netherlands

Tel: +31-43-3881277, Email: s.langie@maastrichtuniversity.nl

39

40 **Abstract**

41 This optimized protocol for the comet assay-based *in vitro* DNA repair assay (including links to
42 instruction videos) is relatively simple, versatile, and inexpensive, allowing the detection of base and
43 nucleotide excision repair activity. Protein extracts from samples are incubated with agarose-
44 embedded substrate nucleoids ('naked' supercoiled DNA), containing specifically induced DNA lesions
45 (e.g., through oxidation, UVC or benzo[a]pyrene-diolepoxide treatment). DNA incisions produced
46 during the incubation reaction are quantified as strand breaks after electrophoresis, reflecting the
47 extract's incision activity. An additional step, supplementing the extract with dNTPs, allows the
48 measurement of ligation activity. Various innovations and optimizations have increased the assay's
49 throughput and enabled the use of various samples (cell models, blood cells, tissues). Once extracts
50 and substrates are prepared, the assay can be completed within two days. This method represents a
51 unique functional measurement of DNA repair activity, with applications in human biomonitoring, *in*
52 *vitro*, *in vivo*, and (clinical) intervention studies.

53

54

55

56

57

58 **Key words:** DNA repair, incision, comet assay, nucleotide excision repair, base excision repair, protein
59 extract, single-cell gel electrophoresis

60

INTRODUCTION

The comet-based *in vitro* DNA repair assay is a modified version of the comet assay (also known as single-cell gel electrophoresis assay) to assess DNA repair activity: a cellular protein extract which contains repair enzymes is incubated with a DNA substrate containing induced lesions, and DNA incisions in the form of DNA strand breaks (SBs) accumulate. It is a relatively simple method for functional measurement of base excision repair (BER) and nucleotide excision repair (NER) activity of different types of samples, with many applications in human biomonitoring, in *in vitro* and *in vivo* studies, as well as in (clinical) intervention studies.

Development of the protocol

The comet assay is a versatile and sensitive method that detects, in its standard version, SBs and alkali-labile sites. The first paper on this single-cell gel electrophoresis assay was published in 1984 by Ostling and Johanson¹. The protocol is simple: briefly, cells embedded in agarose on a microscope slide are lysed to remove membranes and soluble components (including histones) leaving nucleoids (i.e., supercoiled DNA attached at intervals to a nuclear matrix forming loops)². Next, nucleoids undergo alkaline unwinding and electrophoresis. The presence of SBs in the DNA relaxes the supercoiled loops and enables the DNA to migrate towards the anode. The resulting comet-shaped figures, called comet(s), are visualised with a DNA fluorescent dye and fluorescence microscopy. In addition, the enzyme-modified comet assay includes an extra step between lysis and alkaline treatment; i.e. incubation with DNA repair enzymes from bacteria or human cells to gain further information on specific classes of DNA lesions³. For instance, among others, formamidopyrimidine DNA glycosylase (Fpg) detects oxidized purines, formamidopyrimidines (ring-opened adenine or guanine) and ring-opened N7 guanine adducts; human 8-oxo-guanine (8-oxoG) DNA glycosylase (hOGG1) detects oxidized purines and formamidopyrimidines; and T4 endonuclease V (T4endoV) detects dimerised pyrimidines. These enzymes are also used in combination with the comet-based *in vitro* DNA repair assay as ‘incubation reaction controls’.

The comet assay (with and without the inclusion of lesion-specific enzymes) is widely used as a biomarker assay in human population studies and genotoxicity testing (including regulatory toxicology) - primarily to measure DNA damage, but increasingly also to assess the activity of cells for DNA repair. In the original publication, Ostling and Johanson also reported the first experiments to measure DNA repair by simply following the decrease of ionising radiation-induced SBs over time – referred to as a challenge assay or cellular repair assay. However, this approach merely measures the final step in the repair process (i.e., ligation). Still, useful information on the kinetics of NER and BER has been gained by following the removal of pyrimidine dimers or oxidised bases, respectively, using appropriate

enzymes^{4,5}. Then again, this approach is time-consuming and laborious, and therefore not optimal for biomonitoring or intervention studies, which typically require high-throughput processing of many samples.

An alternative *in vitro* approach is based on assessing the ability of repair proteins in a sample extract to recognize and incise substrate DNA that contains induced lesions. The whole-cell extract can be prepared from blood cells, ground tissues or cultured cells, by 'snap-freezing' and subsequent lysis with Triton® X-100. The comet-based *in vitro* DNA repair assay was first devised in 1994 to measure NER and BER activity in a human cell extract⁶. However, over the past two decades, it has been modified and improved, as well as being applied to tissue samples in addition to cell suspensions. **Table 1** gives the comparison between the three main versions of the comet assay, i.e. the standard comet assay, the enzyme-modified comet assay, and the comet-based *in vitro* DNA repair assay.

Variations in the method

The nature of the lesions in the substrate nucleoids defines the repair pathway that is going to be studied. Early on, the comet-based *in vitro* DNA repair assay was applied to extracts from human lymphocytes to measure BER, using cells treated with a photosensitizer (Ro 19-8022) plus visible light to create substrate DNA containing 8-oxoG⁷. The use of the photosensitizer Ro 19-8022 has certain challenges; including that for a long time, it was only available on request from F. Hoffmann-La Roche, and that the irradiation with light may increase the level of DNA SBs. In the past 5 years, potassium bromate (KBrO₃) has been introduced as a novel, alternative DNA damage-inducing agent to prepare substrate cells containing oxidatively damaged DNA lesions. In addition, Fpg has a higher incision rate than hOGG1 on Ro 19-8022 treated cells, whereas these repair enzymes have the same incision rate in KBrO₃-generated substrate cells⁸. So far, KBrO₃-treated substrate cells have been used in only a few studies to assess DNA repair activity in cell cultures⁹, animal tissues¹⁰, and human blood cells¹¹.

The BER-specific *in vitro* DNA repair assay was modified in 2005 for evaluation of NER activity in human lymphocytes, using benzo[a]pyrene-diolepoxide (BPDE) – the active metabolite of the well-studied environmental mutagen benzo[a]pyrene (B[a]P) – to treat substrate cells¹². This type of DNA-damaging agent induces BPDE-DNA adducts which are typically recognized by NER enzymes. The original version of this assay involved the treatment of nucleoids in the gel with BPDE after lysis¹². Recent optimizations, performed by Working Group 5 of the hCOMET COST Action (CA15132), improved the standardization of the assay and successful creation of a batch of BPDE-exposed substrate cells (as explained in Step 1-D). In 2009 the *in vitro* DNA repair assay was further modified by using UVC-treated substrate cells to measure NER¹³. **Box 1** provides an overview of the various modifications of the assay,

including links to the corresponding steps in the protocol. An advantage of the assay is its versatility, illustrated by studies that have adopted alternative DNA damage-inducing agents, such as oxaliplatin, to study repair of DNA cross-links or H₂O₂ and methyl methanesulfonate (MMS) to induce various base modifications^{14,15}.

In most studies, the comet-based *in vitro* DNA repair assays have used protein extracts from cultured or isolated cells (e.g. blood cell fractions) to study DNA repair activity. Various attempts have been made by different laboratories to use the comet-based *in vitro* DNA repair assays with extracts from solid animals tissues, but only a few have succeeded^{16,17}, most being frustrated by low repair activity and/or low detection sensitivity due to the presence of non-specific nuclease activity¹⁸. However, from 2010 onwards, methods for assessing BER and NER activities were developed and optimized for their use with protein extracts from solid tissues. While Langie *et al.* modified both BER and NER assays to measure DNA repair activities from solid tissues of animal origin^{16,18}, human tissue samples were assayed for both repair pathways by Slysikova *et al.*¹⁹. It is our aim to optimize the assays further for use with non-invasively collectable tissues, such as buccal cells and saliva – which are currently used successfully in the standard alkaline comet assay to assess DNA damage levels.

BOX 1: OVERVIEW OF THE VARIOUS MODIFICATIONS OF THE COMET-BASED *IN VITRO* DNA REPAIR ASSAY

Schematic overview of the most used modifications of the comet-based *in vitro* DNA repair assay, including links to the corresponding steps in the protocol.

BER assay

- Substrate cells ➤ [Treatment with photosensitizer Ro 19-8022 + light](#)
- [Treatment with KBrO₃](#)
- Incubation step ➤ [To study incision activity – use extract + buffer B](#)
- [To study synthesis & ligation – use extract + buffer B + dNTPs](#)

NER assay

- Substrate cells ➤ [Treatment with UVC](#)
- [Treatment with BPDE](#)
- Incubation step ➤ [To study incision activity – use extract + buffer N](#)
- [To study synthesis & ligation – use extract + buffer N + dNTPs](#)

Protein extract preparation

- [From frozen cells or PBMC/WBC in freezing medium](#)
- [From fresh or frozen cell pellets](#)
- [From snap frozen tissues](#)

Abbreviations: BER – base excision repair, BPDE – benzo[a]pyrene-diolepoxide, dNTPs – deoxyribonucleotides, KBrO₃ – potassium bromate, NER – nucleotide excision repair, PBMC – peripheral blood mononuclear cells, UVC – ultraviolet (C), WBC – white blood cells.

Principle of the assay

Figure 1 provides a schematic overview of the principle behind the sample extract incubation reaction; it depicts substrate DNA from cells that were pre-treated, as an example, with the photosensitizer Ro 19-8022 plus light for the measurement of BER. In general, a protein lysate is extracted from cells or tissues and incubated with damage-containing substrate DNA, consisting of gel-embedded nucleoids from cells that were pre-treated with DNA damage-inducing agent. Incubation of these substrate nucleoids with cell or tissue extracts allows the initial steps of BER (or NER in the case of UVC- or BPDE-induced DNA lesions) to occur; repair enzymes present in the protein extracts will induce incisions at the site of the DNA lesions in the substrate. These incisions will result in single SBs that can then be determined by the standard alkaline comet assay. Thus, the increased migration of DNA into the tail is proportional to the DNA repair incision activity of the extracts.

This assay essentially assesses the DNA incision activity, measuring the accumulation of DNA SBs, and incision is generally regarded as the rate-limiting step of DNA repair. Therefore, by merely assessing the DNA repair incision activity, it is already possible to study the effect of external and internal factors on an organism's DNA repair activity. In real life, the SBs produced during the incision step are transient, being quickly followed by DNA repair synthesis (long patch synthesis in the case of NER, shorter patches or single nucleotides in BER). *In vitro*, the concentration of deoxyribonucleotides (dNTPs) is too low for this synthesis to occur – as is confirmed by the experimental addition of ATP and dNTPs, which prevented the increase in SBs after UVC irradiation in HeLa cells⁶. Thus, if it is also required to detect DNA synthesis & ligation activity of the extract, a parallel incubation of sample extracts supplemented with dNTPs can be performed (Box 1, and Optional step in parallel to Step 4 presented in detail in Box 5).

Applications of the method

The comet-based *in vitro* DNA repair assay has been used in some cell culture and animal studies – studying the effect of nutrition and ageing – but it is mostly used in human biomonitoring and intervention studies. We previously reviewed the different *in vitro*, *in vivo* animal and human studies where this technique has been applied to measure DNA repair activity²⁰. The text below gives the main messages from this review. In the near future, we also plan to use the assay in genotoxicity testing to unravel the role of DNA repair in the Mode-of-Action (MoA) of potential (non)genotoxic carcinogens. In addition, DNA repair has recently been defined as a key event (KE) in an adverse outcome pathway (AOP) that was submitted to the OECD Extended Advisory Group for Molecular Screening and Toxicogenomics (EAGMST) for internal review²¹ – which may also promote the use of the assay.

Cell culture studies

There are few studies in the literature where the comet-based *in vitro* DNA repair assay has been applied using cell cultures. Most of them studied the (beneficial) effects of nutrients, mainly polyphenols²²⁻²⁴ and other antioxidants^{5,25}, on the DNA repair activity. In a few cases, the effects of therapeutic drugs^{26,27} on DNA repair activity were tested, or the assay was used to unravel underlying disease mechanisms¹⁵.

Animal studies

The first three reports of the use of the comet-based *in vitro* DNA repair assay on animal tissues (rodents and pigs) only came about 8-9 years after the first reports on the assay using human blood samples^{16,17,28}. This slow start was due to the presence of high levels of non-specific incision activity when using protein extracts from tissues, making the measurement of DNA repair in mammalian tissues using the comet-based assay a challenge. The adapted and optimized assay¹⁸ for quantification of BER-associated incision activity in rodent tissues opened opportunities for a wide range of *in vivo* studies, including effects of environmental exposures (such as toxins, dietary factors, and pharmaceutical agents) and physiological processes including growth, development, degenerative diseases, and ageing. The comet-based *in vitro* DNA repair assay has mainly been used to study the effect of ageing or dietary factors in animal tissues^{17,29-31}. However, in recent work by Setayesh *et al.*, the effect of weight-loss strategies on the NER activity in obese mice was studied³².

In 2014, the comet-based *in vitro* DNA repair assay was applied to *Drosophila melanogaster* to measure the DNA repair activity in extracts from different strains, proficient and deficient in DNA repair³³. The *in vitro* approach can provide information about the genetic basis and regulation of specific repair enzymes.

Human studies

Individual DNA repair activity is a valuable biomarker since it has been regarded as a marker of susceptibility to mutation and cancer development. A high repair activity is related to a decrease in the chance of unrepaired damage when cells replicate and so to a decrease in potential mutations.

The comet-based *in vitro* DNA repair assay has been used mainly in human biomonitoring and nutritional intervention studies, but also in occupational and clinical studies³⁴⁻³⁷. In a recent review we give an overview of the use of the comet-based *in vitro* DNA repair assay in various human biomonitoring studies and describe how DNA repair activity can be affected by various external (e.g., chemicals, lifestyle, diet) and internal (e.g., genetics, age, sex) factors³⁸. In occupational studies, while the harmful effect of exposure was clearly recognizable by high levels of various biomarkers of genotoxicity, the effect of exposure on DNA repair activity was not always that straightforward. For

instance, in workers exposed to stone wool, BER activity was unaffected by exposure but was negatively correlated with micronucleus frequency, implying that unrepaired 8-oxoG contributes to micronucleus formation³⁹. However, a study on occupational exposure to asbestos showed that non-exposed women had higher mean BER activity compared with exposed women⁴⁰. In styrene-exposed hand-lamination workers, an exposure-related increase in BER activity and a decrease in SBs was observed, suggesting possible induction of DNA repair enzymes in the course of chronic occupational exposure⁴¹.

In studies investigating DNA repair activity in relation to human diseases, the comet-based *in vitro* DNA repair assay has been used only rarely, mainly on peripheral blood mononuclear cells (PBMC) of study subjects (e.g., patients suffering chronic renal failure⁴², patients with lung⁴³ and colorectal cancer^{44,45}), but also on biopsies from tumour and adjacent non-tumour tissue from colorectal cancer patients^{19,46}. Since the comet-based *in vitro* DNA repair assay to study BER and NER in human solid tissues was optimized only recently^{19,47}, more clinical studies on DNA repair in relation to tissue-specific diseases might be expected to be published in the near future.

Results to date have demonstrated the range of repair activities in a healthy human population – a range far wider than can be explained by genetic polymorphisms. This emphasizes the importance of regulation of repair by environmental and/or intrinsic factors – about which we still know relatively little. Nonetheless, the assay allows the assessment of the intrinsic DNA repair activity, as observed from measurements for the same persons at different time points. The comet-based *in vitro* DNA repair assay is the perfect tool to phenotypically assess the activity of various DNA repair pathways and thereby to further unravel the effect of various modifying factors on the activity as well as investigating the DNA repair activity as an effect modifier in studies on exposures to genotoxic agents.

Comparison with other methods

Inducing DNA damage in cells and monitoring the rate of removal of the lesions over time is the most straightforward approach to measuring DNA repair activity (also known as cellular repair assay). The comet assay, in the standard version to measure DNA damage, has been used with this aim since the very beginning of the assay^{1,48}. Moreover, the use of polymerase inhibitors (e.g., aphidicolin or cytosine arabinoside) increases the sensitivity of the assay⁴⁹⁻⁵¹. Three *in vitro* studies have demonstrated an increase in the rate of removal of oxidized bases or DNA SBs in line with an increase in BER activity estimated by using the comet-based *in vitro* DNA repair assay^{5,23,24}. However, from the logistical point of view, this approach is not very convenient when analysing a large batch of samples.

Several analytical techniques can be used to monitor the removal of the damage over time. A significant positive correlation was observed between the NER comet-based *in vitro* DNA repair assay

and BPDE-DNA adduct removal determined by the ^{32}P -post-labelling assay¹². According to our knowledge, additional comparative studies between analytical techniques (e.g. HPLC or MS) and the comet-based *in vitro* DNA repair assay have not been carried out.

There are other methods to assess DNA repair activity. The unscheduled DNA synthesis (UDS) assay, used for many years, is based on the incorporation of [^3H] thymidine into the DNA after treatment with a genotoxic agent⁵²; it is effective in measuring the repair of UVC-induced damage but less effective in measuring the smaller gaps produced during BER. In any case, this method has been widely criticised, and it is not recommended nowadays.

Different approaches using plasmids have been developed to measure DNA repair activity. In a host cell reactivation (HCR) assay⁵³, a plasmid containing a UVC- or BPDE-damaged reporter gene (e.g., luciferase gene) is introduced into the cells. The activity of the reporter gene gives an estimate of the extent to which the cells have been able to repair the lesion in the plasmid. *In vitro* plasmid-based assays are more common. Plasmids containing DNA lesions are incubated with cell extracts containing repair enzymes. Then, using a standard gel electrophoresis method, nicked plasmids (repaired) can be separated from closed, non-repaired ones⁵⁴. Plasmids can also be incubated with the cell extract in the presence of ^{32}P -labelled dNTPs and the repair measured by their incorporation into the plasmid⁵⁵. In this way, not only incision but the whole repair process is measured. An alternative to the use of plasmids is the use of oligonucleotides constructed with specific lesions and a terminal radioactive or fluorescent tag^{56,57}. To date, no direct comparison of the comet-based *in vitro* DNA repair assay with plasmid or oligonucleotide-based techniques has been carried out. Although such methods have been applied in human biomonitoring studies, especially by Paz-Elizur *et al.*⁵⁸⁻⁶⁰ and Leitner-Dagan *et al.*^{61,62}, the number of studies in which these techniques were applied is limited. The comet assay, on the other hand, has been used as an *in vitro* DNA repair assay more often⁶³.

DNA repair has also been measured in terms of the level of transcription of DNA repair-related genes. However, an important limitation of this approach is that post-transcriptional regulation and epigenetic changes are not taken into account. In fact, there seems to be little correlation between transcription levels and repair enzyme activities^{60,64-66}. Thus, DNA repair phenotyping is a more direct measurement than genomic, transcriptomic and proteomic approaches.

The *in vitro* comet-based approach to measuring DNA repair activity became increasingly popular due to several positive aspects. Above all, it opens up the possibility to study DNA repair activity in diverse biological material, unlike cellular repair assays in which the removal of damage over time is monitored in cells in culture. The cell extract to be used in the *in vitro* assay can be prepared from virtually any

tissue. Moreover, frozen materials (cells in freezing medium, frozen cell pellets, or tissues) can be used to prepare protein extracts, which make it logistically a more attractive assay.

As compared to other DNA repair assays such as UDS, HCR, and plasmid- or oligonucleotide-based incision assays, the comet-based *in vitro* DNA repair assay detects the effect in nucleoids (condensed/supercoiled DNA, as occurs in cells), and it is not necessary to use radioactively-labelled material. Moreover, the density of lesions in the nucleoid is low, which may represent a theoretical advantage over other assay designs since it more closely resembles the environment that the repair enzymes encounter *in vivo*. However, it is still an artificial environment; ‘naked’ supercoiled DNA is not the natural substrate for repair.

Furthermore, the comet-based *in vitro* DNA repair assay is far less laborious and time-consuming compared to the cellular repair assay and can, therefore, be performed on many samples in parallel on a large scale. The practical advantages of the comet-based *in vitro* DNA repair assay are low-cost, simplicity, and versatility. It involves one simple incubation step, and the results on DNA incision activity are obtained within two days. It is therefore well suited for biomonitoring or intervention studies, or for the screening of new chemicals and therapeutics, which typically require high-throughput processing of many samples.

Experimental design

Overview of the assay

The whole procedure of the comet-based *in vitro* DNA repair assay can be divided into six major steps, as described in **Figure 2**. The most crucial steps are also demonstrated in the associated films (<https://www.youtube.com/playlist?list=PLEVxCdaQpbj1GDqGUHgWiaBy9eVTUZOzX>). Steps 1 and 2 have to be performed on day 0, Steps 3-5 on day 1, and Step 6 on day 2.

Step 1: Various types of cells can be used for preparing the substrate DNA – cultured cells or freshly isolated PBMC are convenient. The aim is to produce non-exposed cells with a low background level of DNA damage and exposed cells with a sufficiently high level of specific DNA damage for the enzymes in the extract to work on (obeying the biochemical principle that the lesions in the substrate DNA should be present in excess). (Video instructions: <https://youtu.be/awtdmFBI1WA>)

Checkpoint 1: DNA repair activity may vary between experiments for various reasons, including the amount of DNA lesions induced in the substrate cells. The ratio of specific DNA damage (8-oxoG or UVC-photoproducts/BPDE products) to non-specific damage (SBs and alkali-labile sites) should be verified in preliminary experiments. Such verifications involve the incubation of an endonuclease

specific for the type of DNA damage that is introduced in the substrate cells (Table 3). For expected results and recommended levels of DNA damage in both non-exposed and exposed cells, see Figure 9 (for Ro19-8022 and KBrO₃ exposure), Figure 10 (for UVC exposure), and Figure 11 (for BPDE exposure).

Step 2: To prepare protein extracts, various starting materials can be used: PBMC, cultured cells, animal and human tissues. It is advisable to prepare extracts from all the samples at roughly the same time, or at least in large batches, to reduce the risk of batch variations affecting results. (Video instructions: <https://youtu.be/VHRHwkJFIDw>)

Check point 2: Protein concentration should ideally be measured prior to the reaction so that all the extracts can be diluted to and used at the same concentration on the day of experiment. Retrospective normalization of the activity according to protein concentration (when extracts are used at different concentrations in experiments and results are adjusted for protein concentration afterwards) is not recommended because protein concentration and activity as measured in the assay are not proportionally related⁶⁷.

Step 3: When embedding the cells in agarose gel, the final concentration is essential as the migration of DNA and sensitivity of the assay depend on the density of the gel. There are several procedures for performing lysis of the cells in the comet assay. For the comet-based *in vitro* DNA repair assay, a lysis time of 1h is recommended. (Video instructions: <https://youtu.be/T42JOvD2MnE>)

Step 4: At this point in the protocol, substrate agarose-embedded nucleoids, both non-exposed and exposed, are going to be incubated with either reaction buffers or sample extracts (containing DNA repair enzymes). During the incubation, DNA repair enzymes contained in the sample extract induce DNA SBs at the sites of specific DNA lesions in the substrate nucleoids (8-oxoG for BER, or UVC- or BPDE-induced lesions for NER) (Figure 1).

This is one of the most critical steps of the assay and standardization is necessary regarding the time of incubation and concentration of the extracts (see Material setup). Several experimental controls should be included in the assay for the correct interpretation of the results (described below in section “Controls”). (Video instructions: <https://youtu.be/GzghrROzD64>)

Step 5: As a result of the incisions (i.e., SBs), the DNA will be drawn towards the anode forming a comet-like image. The proportion of total DNA in the comet tail reflects the DNA repair activity of the sample extract, which means more DNA incision activity will result in more DNA in the comet tail. The following steps comprise the neutralization and washing of the microscope slides. (Video instructions: <https://youtu.be/kvgZ7O25kXo>)

Step 6: For visualization of the comets, various dyes can be used. The use of tail intensity (TI, % tail DNA) is advised to express the results. However, other primary comet assay descriptors (e.g. tail moment or visual score) can be used to calculate the final DNA repair incision activity.

Controls

Positive and negative controls

It is essential to document the reliability of the comet-based *in vitro* DNA repair assay by analysing control samples in the validation process and on-going experiments. However, there is a lack of experimental controls - chemical or physical exposures - that consistently have been shown to alter the DNA repair activity without causing cytotoxicity or cell death. Instead, it is possible to use repair-deficient cells or tissue samples as negative controls. Ogg1 knockout fibroblasts and mouse tissues are useful sources of repair-deficient extract in the repair assay on Ro 19-8022- or KBrO₃-exposed substrate cells^{7,68}. Fibroblast cell lines from xeroderma pigmentosum complementation group A and C can be used for the repair assay using UVC and BPDE exposure¹². Heat inactivation of repair extracts is a simpler solution to generate a negative control if DNA repair-deficient cells or tissues cannot be obtained^{17,18,69}. There is currently no “true” positive control in the sense that certain cells have higher than normal DNA repair activity. The development of knock-in cells is theoretically possible for at least some repair pathways, but it has not been common practice to do so.

Internal experimental controls

The comet-based *in vitro* DNA repair assay uses internal experimental controls, which are also used in the calculation of the repair-related DNA incision activity or simply to assess if the assay was performed well. These controls assess the incisions/cleavage in nucleoids from non-exposed or exposed substrate cells, incubated with incubation reaction buffer or sample protein extract:

i) the “**background control**” is non-exposed substrate cells incubated with the incubation reaction buffer to check the basal level of DNA damage in the substrate DNA (Figure 3 – yellow; or Figure 2, Step 4, Microscope slide 1);

ii) the “**treatment control**” is exposed cells incubated with the incubation reaction buffer to reveal the presence of non-specific DNA SBs or alkali-labile sites resulting from the exposure with the damaging agent (Figure 3 – green; or Figure 2, Step 4, Microscope slide 2);

iii) the “**specificity control**” is non-exposed substrate cells incubated with the sample’s protein extract to check for non-specific incision or cleavage activity (Figure 3 – blue; or Figure 2, Step 4, Microscope slide 3);

iv) the “**incubation reaction control**” is exposed substrate cells incubated with a lesion-specific enzyme (Figure 3 – red; or Figure 2, used in Checkpoint 1 and Step 4).

Concerning the latter (iv), it has been common practice to consider the Fpg- or hOGG1- treatment as incubation reaction controls for the KBrO₃ and Ro 19-8022 + light exposed substrate cells^{28,47,70}. T4 endonuclease V has been used as an incubation reaction control for UVC-irradiated substrate cells^{13,47}. However, there are currently no enzymes or crude extracts available that can be used as incubation reaction control for BPDE-generated substrate cells or any other type of bulky DNA adducts that are used to assess NER activity.

Assay setup

Figure 3 illustrates an example of a potential assay setup in which protein extracts of 3 different samples are assessed for their BER and NER incision activity. When preparing the required number of slides, it is important to keep in mind to include the assay controls that were described above. In the scheme, the gels are randomized (e.g. the gels with UVC-exposed substrate cells are in different places on duplicate slides). Alternatively, one can simply put the duplicate gels on the same slide. An example of a setup when using a higher throughput 12-gel system can be found in Box 2.

Limitations

It is worth emphasizing that the repair pathway studied is defined by the kind of damage introduced in the substrate DNA, but there is ‘cross-talk’ between pathways. For instance, in the case of substrate cells containing KBrO₃-induced DNA lesions, it is not absolutely clear if the assay measures the overall BER activity or just the repair incision activity at oxidatively damaged DNA. Nonetheless, the various assays have been optimized by using knock-out cells or tissues for either BER or NER genes, confirming the specificity of the assays^{7,12,18}.

To prevent unreliable results, the presence of haemoglobin and bilirubin during extract preparation should be avoided as they interfere with quantification of the protein concentration, and so may lead to overestimation of the protein concentration of the extract^{71,72}. Therefore, it is advisable to use a protein assay, such as the Lowry-based BIO-RAD DC Protein Assay Kit, which measures protein concentrations at 650–750 nm. At these wavelengths, the absorbance of haemoglobin (high absorption at ~250–600 nm) and bilirubin (high absorption at ~400–500 nm) is negligible, especially when samples are well diluted^{18,73,74}.

If the protein concentration of the cell/tissue extract and the incubation time are not optimized for each specific cell type, it is possible that only low repair rates are detected. Sometimes non-specific

incisions can even exceed the specific incisions produced by repair enzymes, leading to negative repair rates. Therefore, it is important to optimize both the time of incubation and the protein concentration of the extract (see Equipment setup and Checkpoint 2).

The limitation of any comet-based assay relates to the assessment of large numbers of samples. Traditional practice involves processing of comets in relatively large gels on the microscope slides (with one or two gels per slide), which limits the number of samples that can be run in one experiment. This is ameliorated by the development of high throughput versions of the assay, where 12 mini-gels, instead of 2, are run on one microscope slide (Box 2). However, even these gels require manual scoring of the comets in the microscope, which sets a limit to the number of samples that can be processed per working day. Automatic systems have been developed for the identification and scoring of comets, but few researchers appear to use them.

Individual DNA repair activity is regarded as a marker of susceptibility to genotoxic agents, on the assumption that a high intrinsic repair activity will be protective. However, it is still not clear whether a high repair activity might be induced in response to, and therefore indicative of, exposure to DNA damaging agents⁷⁵ - though this is a general limitation for all *in vitro* DNA repair assays. In any case, it is undoubtedly of value to gather information about individual repair activity alongside DNA damage measurements, since the two are intimately connected. The steady-state level of DNA lesions in a cell is determined by the damage input and the capacity of the cell to repair the damage.

MATERIALS

Biological materials

- Cell cultures – using cultured cells (in monolayer or suspension culture) is the most straightforward way to create substrate cells. In addition, cultured cells can be used to prepare cellular protein extracts.
- Animal samples – blood (WBC, PBMC) and different tissues can be used for preparing protein extracts.
- Human samples – blood (WBC, PBMC) can be used to create substrate cells, while blood and different tissues (i.e. biopsies, potentially also buccal cells and saliva) can be used to prepare protein extracts.
- *Drosophila melanogaster* - *Drosophila* larvae cells (neuroblasts, haemocytes and anterior midgut cells) – can be used to prepare sample extracts (for details, see Box 3)

Δ CRITICAL Various cell types, used to create substrate cells, can show different levels of background DNA damage. Similarly, DNA repair activity and non-specific nuclease activity vary with cell or tissue type. Therefore, before each set of experiments with particular cells, it is essential to check the background levels of DNA damage in the substrate cells and titrate the protein concentration of the sample extracts (for details, see the Checkpoint 1 and 2).

Reagents

For all the reagents mentioned below the most commonly used provider is mentioned, though reagents purchased from other providers should perform equally well.

General reagents

- Agarose, normal melting point (NMP) (Sigma-Aldrich, cat. no. A4718)
- Agarose, low melting point (LMP) (Sigma-Aldrich, cat. no. A9414)
- Phosphate buffered saline (PBS) (Sigma-Aldrich, cat. no. P4417)
- Triton® X-100 (Sigma-Aldrich, cat. no. X100)
- 4-(2-Hydroxyethyl) piperazine-1-ethanesulfonic acid (HEPES) (Sigma-Aldrich, cat. no. H3375)
- Ethylenediaminetetraacetic acid disodium salt dihydrate (EDTA- $\text{Na}_2 \cdot 2\text{H}_2\text{O}$) (Sigma-Aldrich, cat. no. E5134)
- DL-Dithiothreitol (DTT) (Sigma-Aldrich, cat. no. D9163)
- Glycerol (Sigma-Aldrich, cat. no. G5516)
- Trizma® base (Sigma-Aldrich, cat. no. T1503)
- Potassium chloride (KCl) (Sigma-Aldrich, cat. no. P3911)

- 452 • Sodium chloride (NaCl) (Sigma-Aldrich, cat. no. S9888)
- 453 • Potassium hydrochloride (KOH) (Sigma-Aldrich, cat. no. P5958)
- 454 • Sodium hydroxide (NaOH) (Sigma-Aldrich, cat. no. 795429) **! CAUTION** NaOH is caustic
- 455 • Bovine Serum Albumin (BSA) (Sigma-Aldrich, cat. no. A2153)
- 456 • Ethanol 96% (Merck Millipore, cat. no. 159010)
- 457 • Liquid nitrogen
- 458 • Adenosine 5'-triphosphate (ATP) (Sigma-Aldrich, cat. no. A1852)
- 459 • Deoxyribonucleotides (dNTPs) (ThermoFisher, cat. no. R0181)
- 460 • The Lowry-based BIO-RAD DC Protein Assay Kit using bovine serum albumin as a standard and
- 461 controlling for the presence of Triton® X-100, DTT, and EDTA (BioRad, cat. no. 500-0116)

462 ***Reagents for cultivation, freezing and counting the substrate cells***

- 463 • Cell culture medium (depending on the cells used)
- 464 • Dimethyl sulfoxide (DMSO) (Sigma-Aldrich, cat. no. 41639)
- 465 • Foetal Bovine Serum (FBS) (Sigma-Aldrich, cat. no. TMS-016)
- 466 • Other; depending on the cells used, cell culture medium may need the addition of some
- 467 complements (e.g., non-essential amino acids, glutamine, penicillin/streptavidin).

468 ***Reagents for preparation and checking of substrate cells***

- 469 • Photosensitiser Ro 19-8022 - for preparing BER substrate (CAS 104604-66-2, can be obtained
- 470 from e.g. Chiron – cat. no. C8504.19-1-DS, Pharmaffiliates – cat. no. PA 27 00232) **! CAUTION**
- 471 Genotoxic, wear protective gloves.
- 472 • Potassium bromate (KBrO₃) - for preparing BER substrate (Merck, cat. no. 104912) **! CAUTION**
- 473 Carcinogenic, toxic, wear protective gloves.
- 474 • Benzo(a)pyrene-r-7,t-8-dihydrodiol-t-9,10-epoxide(+/-) (anti) (BPDE) - for preparing NER
- 475 substrate (Bio-connect BV., cat. no. MBS6101688). **! CAUTION** Carcinogenic, mutagenic, wear
- 476 protective gloves.
- 477 • Formamidopyrimidine DNA glycosylase (Fpg) - Incubation reaction control for BER (New
- 478 England Biolabs, cat. no. M0240S, or obtained from NorGenoTech)
- 479 • Human 8-oxoguanine DNA glycosylase (hOGG1) – alternative incubation reaction control for
- 480 BER (Trevigen, cat. no. 4130-100-EB; or obtained from NorGenoTech)
- 481 • T4 Endonuclease V (T4endoV) - Incubation reaction control for NER (New England Biolabs, cat.
- 482 no. M0308S)

Reagents for comet visualization

Several DNA fluorescence dyes are suitable, however, the most commonly used are:

- SYBR® Gold (ThermoFisher, cat. no. S11494) ! CAUTION Potential mutagen, wear protective gloves.
- SYBR® Green (ThermoFisher, cat. no. S7567) ! CAUTION Potential mutagen, wear protective gloves.
- Ethidium bromide (ThermoFisher, cat. no. 17898) ! CAUTION Mutagenic activity, wear protective gloves.
- DAPI (ThermoFisher, cat. no. D1306) ! CAUTION Mutagenic activity, wear protective gloves.

Other newly developed dyes, such as GelRed®, can be used as well.

Equipment

Common equipment and consumables to perform cell culture or to collect human/animal samples are needed. Moreover, general laboratory equipment and consumables are required (e.g., microwave oven, freezers, fridge, pH meter, cooled centrifuge, plastic tubes, vortex, plastic tips, pipettors, micropipettes). Special equipment and consumables needed for the comet-based *in vitro* DNA repair assay can be obtained from various providers. Although certain providers are recommended, the protocol works with most equipment.

- Microscope slides - standard microscope slides with frosted end are used
- GelBond® films (Lonza, cat. no. 53734) can be purchased for mounting of the gels instead of using microscope slides
- 20x20 mm coverslips, or 22x22 mm coverslips to mould gels
- 24x60 mm coverslips
- 500W tungsten halogen lamp - for activation of the photosensitiser Ro 19-8022
- Germicidal UVC lamp - for induction of UVC-induced damage in substrate DNA
- UVC-dosimeter
- Mr Frosty® (Nalgene, VWR cat. no. 479-3200) freezing container, or tick walled (min. 1 cm) polystyrene box – to slowly freeze substrate cells
- Pestle and mortar
- Hammer
- Nanodrop or plate reader – to quantify protein concentrations
- Microtube pestles - for homogenisation of tissues
- Water bath or thermoblock

- Staining (Coplin) jars - for cell lysis and slide washing
- 12-Gel Comet Assay Unit (NorGenoTech)
- Metal trays or plates – to keep slides cold and prevent enzyme reaction to start
- Incubator + moist box - for extract-substrate incubation (alternative is a heating plate or 'slide moat' purchased from Boekel Scientific)
- Large-bed horizontal gel electrophoresis chamber
- Power supply (one that reaches 1-2 Amp is advised, e.g. obtained from Consort)
- Epifluorescence microscope and filter set for green-light excitation, Charge-coupled device (CCD) camera (8-bit black-and-white camera is adequate); high sensitivity and high pixel density are preferred
- Optional: i) peristaltic pump to recirculate the electrophoresis solution (e.g., there are cheap peristaltic pumps made for aquariums); and ii) recirculating chiller to cool the platform of the electrophoresis tank.

Software

- For scoring comets, computer-assisted image analysis is recommended using commercially available software which gives the most reproducible results. Examples of scoring software: Comet assay IV (Instem), Comet Analysis software (Trevigen), Lucia Comet Assay™ software (Laboratory Imaging), Metafer (Metasystem).
- Several scoring programs are freely available, among which Casplab and CometScore showed a good agreement with the Comet assay IV Software (Instem), while OpenComet (plugin of ImageJ) showed the least agreement – especially when only samples with %Tail DNA<15% were analysed [unpublished data generated by Working Group 5 within the hCOMET COST Action CA15132].
- Alternatively, visual scoring – classifying comets into 5 classes based on the amount of DNA in the tail⁷⁶ – has shown good agreement with commercially offered software⁷⁷, but it needs proper training to classify the comets correctly and objectively.

Reagents setup

Solutions

- Cell freezing medium (for freezing cells): DMEM, 10% (vol/vol) foetal bovine serum (FBS), 10% (vol/vol) DMSO. Mix 8 mL of DMEM, with 1 mL foetal bovine serum and 1 mL DMSO. Prepare fresh on the day of use. If needed it can be stored at 4°C up to 24h. *Note:* The proportion of FBS in the freezing medium will depend on the cell type used.

- 1% (wt/vol) NMP agarose (for pre-coating slides): Dissolve 1 g NMP agarose in 100 mL distilled water (or proportional volume), microwave to dissolve, and cool to about 50-60°C in a water bath. One hundred millilitres are sufficient to coat about 75-100 microscope slides. 1% NMP agarose should always be made up fresh.
 - 0.7% (wt/vol) LMP agarose in PBS (for embedding cells): Dissolve 0.35 g of LMP agarose in 50 mL PBS, microwave to dissolve, make aliquots of 2-5 mL and store at 4°C. Before use, microwave or submerge the aliquot in boiled water to melt the agarose and then cool to 37°C (in water bath or thermoblock).
- Δ CRITICAL** It is best not to reheat LMP agarose aliquots (otherwise evaporation can cause a significant increase in concentration).
- Δ CRITICAL** The final concentration of the LMP agarose gel, after mixing with the substrate cells, should be $\leq 0.8\%$ (wt/vol), because higher concentrations reduce the sensitivity of the assay.
- Buffer A (extraction buffer): 45 mM HEPES, 0.4 M KCl, 1 mM EDTA- Na_2 , 0.1 mM DTT, 10% (vol/vol) glycerol. Prepare 100 mL. Dissolve 1.07 g HEPES, 2.98 g KCl, 37.20 mg EDTA- $\text{Na}_2 \cdot 2\text{H}_2\text{O}$, 1.54 mg DTT into 90 mL of distilled water. Add 10 mL of glycerol. Adjust to pH 7.8 with 10 M KOH (dissolve 280.55 g in 0.5 L distilled water). Store frozen (-20°C) as 1-2 mL aliquots. Stable for at least 6 months.
 - Buffer A/1% (vol/vol) Triton® X-100: Prepare 1% Triton® X-100 in buffer A: add 10 μL of Triton® X-100 to 990 μL of buffer A. Store frozen (-20°C) in 1 mL aliquots (for use in single experiment). Stable for at least 6 months.
 - Buffer A/0.25% (vol/vol) Triton® X-100 (for background control incubation): Prepare 0.25% Triton® X-100 in extraction buffer A: add 2.5 μL of Triton® X-100 to 997.5 μL of buffer A. Store frozen (-20°C) in 0.5 mL aliquots (for use in single experiment). Stable for at least 6 months.
 - Lysis solution: 2.5 M NaCl, 0.1 M EDTA- Na_2 , 10.0 mM Trizma® base: Dissolve 146.10 g NaCl, 37.22 g EDTA- $\text{Na}_2 \cdot 2\text{H}_2\text{O}$, 1.21 g Trizma® base into 1 L of distilled water. Adjusted to pH 10 with 10 M NaOH (dissolve 200 g of NaOH in 0.5 L distilled water). Prepare 1 L. Will be stable for at least 6 months when stored at 4°C. Before use, add 1 mL of Triton® X-100 per 100 mL.
 - Buffer B (washing buffer after lysis and incubation reaction buffer for BER): 40 mM HEPES, 0.5 mM EDTA- Na_2 , 0.2 mg/ mL BSA, 0.1 M KCl: Dissolve 9.53 g HEPES, 7.45 g KCL, 0.19 g EDTA- $\text{Na}_2 \cdot 2\text{H}_2\text{O}$, 0.2 g BSA in 1 L distilled water. Adjusted to pH 8 with 10M KOH (dissolve 280.55 g in 0.5 L distilled water). We advise to prepare 500 mL of 10x concentrated stock and freeze (-20°C) in 50 mL tubes (to use for washing slides after lysis) and in 1 mL aliquots (to use as incubation reaction buffer). Stable for at least 6 months. Dilute 10x in distilled water on the

day of use. *Note:* The diluted buffer B could be stored at 4°C for use in a second assay within the same week.

- Buffer N (washing buffer after lysis and incubation reaction buffer for NER): 45 mM HEPES, 0.25 mM EDTA-Na₂, 0.3 mg/ mL BSA, 2% (vol/vol) glycerol. Dissolve 10.72g HEPES, 0.093 g EDTA-Na₂·2H₂O, 0.3 g BSA into 980 mL distilled water. Add 20 mL of glycerol. Adjusted to pH 7.8 with 10 M KOH (dissolve 280.55 g in 0.5 L distilled water). We advise to prepare 500 mL of 10x concentrated stock and freeze (-20°C) in 50 mL tubes (to use for washing slides after lysis) and in 1 mL aliquots (to use as incubation reaction buffer). Will be stable for at least up to 6 months. Dilute 10x in distilled water on the day of use. *Note:* The diluted buffer N could be stored at 4°C for usage in a second assay within the same week.
- Electrophoresis solution: 0.3 M NaOH, 1 mM EDTA-Na₂: Mix 60 mL 10M NaOH (dissolve 200g of NaOH in 0.5 L of distilled water) and 10 mL 200mM EDTA-Na₂ (dissolve 74.45 g of EDTA-Na₂·2H₂O in 1 L distilled water) in 1930 mL of cold distilled water. Store at 4°C for up to one week.
- Neutralising solution: 1xPBS. Store at 4°C or according to manufacturer instructions.
- TE buffer (for SYBR® Gold and SYBR® Green) – Mix 10 mL of 1 M Trizma® base (dissolve 60.57 g in 0.5 L distilled water) and 2 mL of 0.5 M EDTA-Na₂ (dissolve 18.61 g EDTA-Na₂·2H₂O in 100 mL distilled water) in 988 mL of distilled water. Prepare 1L and store at room temperature (approx. 22°C). Will be stable for at least up to 6 months. Alternatively, it is possible to use TBE or TAE buffer as recommended by a provider.

Materials setup

Pre-coating microscope slides

- 1) Prepare 1% (see Reagents setup) NMP agarose solution in H₂O in microwave and keep it at 50-60°C in water bath.
- △ **CRITICAL** To prevent boiling over, you can put the microwave at lowest power for longer time.
- 2) Dip the slides into the gel until the frosted part.
- 3) Wipe one side of the dipped slide and put the slide flat to dry on a heating plate or overnight on the bench. Remember to indicate with a mark on the frosted part, which side of the slide is the coated one.
- 4) Store them in boxes at room temperature. They can be used for at least up to 12 months.

? TROUBLESHOOTING.

Equipment setup

Most of the equipment does not require any special setup, apart from those mentioned below. For your information, particular setups are also demonstrated in the associated video protocol (video <https://www.youtube.com/playlist?list=PLEVxCdaQpbj1GDqGUHqWiaBy9eVTUZ0zX>).

• Equipment setup for exposure of substrate cells

The exposure of substrate cells to either Ro 19-8022 or UVC requires some specific instructions for the setup (see **Figure 4**, and <https://youtu.be/awtdmFBI1WA>).

1) To perform the exposure to Ro 19-8022 + light, a 500 Watt lamp needs to be mounted on a stand about 33cm above the cells on ice to expose the cells for 5min.

Note: Alternatively, a 2000 Watt lamp at 33cm from the cells can be used for 2 min.

2) For the UVC exposure, you can use any UVC lamp (even those in a PCR hood or cell culture cabinet).

(A) First measure the intensity of the lamp in mW/cm² with UVC Radiometer.

(B) The time of exposure needed to achieve 1-2 J/m² can be calculated using the next formula:

Time (seconds) = (E (mJ/cm²))/(I (mW/cm²)); with (I) for the intensity measured by the UVC Radiometer and the energy (E) is recommended to 0.1-0.2 mJ/cm².

(C) If the measured dose is too high, prepare a box or other device with layers of gauze to reduce the intensity until you achieve a measurable timing. E.g. in our hands 6 layers of gauze gave an intensity of 0.0040 mW/cm², leading to an exposure time of 25 seconds.

• Electrophoresis setup

Since the duration of electrophoresis and the electric potential (voltage drop across the electrophoresis tank platform) are the most important drivers of DNA migration⁷⁸, these should be measured and standardized for all your experiments. (Video instructions: <https://youtu.be/kvgZ7O25kXo>)

1) Measure the width of the platform in the electrophoresis tank.

2) Add sufficient volume of electrophoresis solution to cover the microscope slides with at least 5 mm of liquid.

3) Switch on the power supply and measure the voltage over the platform (holding the electrodes on either end of the platform).

△ CRITICAL Ensure that a power supply is used which can supply the output current at the constant voltage and with sufficient volume of liquid (a power supply that can reach 1-2 Amperes should do the job).

4) Calculate the Electric Potential * Time (EPT) value (dimension: V/cm*min) for your setup. *Note:* for comet-based *in vitro* DNA repair assays we recommend adopting a V/cm*min=30.

5) This EPT value (in this case 30 V/cm*min) should be used for all your experiments and reported in publications, to allow comparison between experiments and laboratories.

Optional:

- Use an external peristaltic pump to recirculate the electrophoresis solution. The advantages are: (1) stable conditions allowing more precise measurement of the electric potential; (2) more stable temperature during electrophoresis; and (3) (probably) reduced variations in the local electric potential.
- Use a recirculating chiller to cool the platform of the electrophoresis tank. Alternatively, the electrophoresis tank can be put in a cold room, dedicated fridge or even put on ice.

Before starting the experiment, it is essential to optimize/titrate the protein concentration of the extract as well as the lesion-specific enzymes, and to determine a suitable incubation time.

Optimization of incubation time

The optimal time of incubation should also be established; the incubation time is normally around 10-30 minutes. The preferred incubation time for the extracts should be the one that allows the detection of the DNA repair incision activity in the linear part of the activity-with-time curve (Figure 5). A typical curve shows an initial linear increase in DNA incision activity after which it reaches a plateau. Ideally an incubation time would be selected that is still on the linear part of the time-incisions curve, detecting a high enough number of incisions before reaching the plateau.

Titration of the lesion-specific enzyme

As indicated above, lesion-specific enzymes are used as “incubation reaction controls”, but also to check the levels of induced DNA lesions in the substrate cells (see Checkpoint 1). However, it has been shown that the enzymes from different producers differ both in their activity and specificity towards nucleobase lesions⁸. Therefore, it is important to perform titration curve experiments to achieve optimal conditions for the enzyme treatment. The optimal concentration elucidated from the titration experiments should detect the maximum enzyme-sensitive sites without inducing non-specific SBs. To do so, substrate nucleoids containing the correspondent lesions and substrate without lesions should

be used (as explained in Checkpoint 1). Muruzabal *et al.*⁷⁹ have described how to perform and interpret titration experiments.

Optimization/titration of the protein concentration of the extract

In addition to the incubation temperature, the incision activity of an extract is also dependent on its protein concentration. Although in the case of extracts from cells (i.e. either cell lines or PBMC/WBC) assay conditions are set according to cell concentration, cell counts are not always reliable (and some cells are invariably lost during centrifugation etc.). Therefore, we recommend measuring the concentration of protein in each extract. In the case of extracts from tissue, the protein estimation is essential. To identify an “optimal” protein concentration giving maximal discrimination between lesion-specific incision and non-specific nuclease action of the extract, we recommend running a titration experiment with different protein concentrations of extracts isolated from the cell type/tissue of interest, before starting the main experiments. In **Figure 6**, the general concept of selecting the “optimal” protein concentration (Figure 6 (B)) is shown; too low protein concentrations will yield low background, but also low repair specific incisions (Figure 6 (A)). At high concentrations, on the other hand, non-specific nucleases may increase the background (Figure 6 (C)) leading to suboptimal analysis of the repair rates. Although each lab should optimize their conditions, guidance can be given to the approximate concentration of protein needed as shown in Table 2.

Note: An option to reduce the non-specific activity of the extracts involves adding aphidicolin to the extract from a stock in DMSO, to a final concentration of 1.5 μ M. Aphidicolin (a DNA polymerase inhibitor), when added to the protein extract has been shown to block non-specific nucleases in the BER assay¹⁸.

PROCEDURE

Step 1: Preparation of substrate cells – day 0, • **Timing** 2-6h (depending on the DNA damaging agent)

Exposure of substrate cells to DNA-damaging agents

Note: any cell type can be used, but it is advisable to use cells in suspension to avoid trypsinization and centrifugation steps (see Box 4 for advice on cell types).

1) Prepare the desired number of cells. Prepare enough flasks/dishes for both exposed and non-exposed cells (to serve as treatment controls).

(A) Prepare cell suspension in cell culture medium without FBS.

(i) PBMC are obtained from venous blood and isolated using a standard density gradient centrifugation method. Cell lines that grow in suspension can also be used.

(ii) Count a sample of the cell suspension.

(iii) Centrifuge cells at about 150-300xg, for 5 min.

(iv) Wash cells with PBS and spin again.

(B) Prepare adherent cell flasks.

(i) Cells are grown in a flask or dish in culture medium to near confluence.

(ii) Before exposure, remove medium and wash cells with PBS.

2) Treatment of the cells with DNA damage-inducing agent

(A) *Ro 19-8022 exposure - Induction of oxidative lesions to study BER*

(i) Resuspend the pellet with cold PBS containing photosensitiser Ro 19-8022 or add it to the flask with adherent cells. Generally, the final concentration of 1-2 μ M.

▲ CRITICAL Avoid excessive light by wrapping the tube containing Ro 19-8022 solution in aluminium foil.

! CAUTION Ro 19-8022 is a carcinogenic agent. Wear protective gloves.

(i) Place cells on ice, 33 cm from a 500 W tungsten halogen lamp and irradiate for 5 min. *Note:* plastic is transparent to this visible light.

(ii) Remove Ro 19-8022 solution and wash cells as described below.

(iii) Also prepare control cells, with no photosensitizer in PBS, exposed to light only (non-exposed substrate cells).

(B) *KBrO₃ exposure - Induction of oxidative lesions to study BER*

(i) Resuspend the cells in culture medium and keep at 37°C.

- 733 (ii) Prepare a 10x stock solution of the final KBrO_3 concentration (e.g. 50 mM but may depend
734 on the cell type) in 37°C warm cell culture medium (without FBS). *Note*: the stock solution is
735 only ten-times lower than the water solubility limit of KBrO_3 .
736 **! CAUTION** KBrO_3 is a carcinogenic agent. Wear protective gloves.
- 737 (iii) Mix the cell suspension with the KBrO_3 stock suspension in with 9:1 ratio (e.g. the final
738 concentration of KBrO_3 in the cell suspension will be 5 mM).
739 **Δ CRITICAL** DNA damage induced by KBrO_3 is largely dependent on intracellular GSH.
740 Therefore, it is advised always to run a dose-response experiment first.
- 741 (iv) Incubate the cells for 1 h at 37°C.
- 742 (v) Remove KBrO_3 solution and wash cells as described below.
- 743 (vi) Also prepare control, non-exposed cells that have been incubated only with cell culture
744 medium.

745 *(C) UVC exposure – induction of thymidine dimers to study NER*

- 746 (i) Start with cells in cold PBS.
- 747 (ii) Place the dish under a measured UV source and irradiate with 1-2 J/m² of UVC.
- 748 **Δ CRITICAL** Irradiate without a dish lid, as plastic reduces UVC exposure.
- 749 **? TROUBLESHOOTING**
- 750 (iii) Also prepare control, non-exposed cells.

751 *(D) BPDE exposure – induction of bulky BPDE-DNA adducts to study NER*

- 752 (ii) Resuspend the pellet in cold PBS containing BPDE or add it to the flask with adherent cells.
753 Generally, the final concentration is 1-3 μM.
- 754 **Δ CRITICAL** Substrate cells should be prepared with the active metabolite BPDE and not
755 with the parent compound B[a]P, because during the metabolism of B[a]P reactive oxygen
756 species may be formed leading to DNA damage that is typically repaired by BER enzymes
757 instead of NER.
- 758 **! CAUTION** BPDE is a carcinogenic agent. Wear protective gloves.
- 759 (iii) In parallel non-exposed cells will be treated with vehicle control DMSO.
- 760 **Δ CRITICAL** The DMSO concentration should be kept as low as possible; <0.5% is advised to
761 avoid toxic effects and to ensure low background DNA damage.
- 762 (iv) Incubate the cells for 30 min at 4°C.
- 763 (v) Remove BPDE or DMSO solution and wash cells as described below.
- 764

3) After treatment, adherent cells need to be washed with cold PBS and trypsinized. In case of cells in suspension go directly to the next step.

△ CRITICAL Avoid long trypsin treatment since this can increase background damage.

4) Transfer the cells to tubes and centrifuge for 5 min at 4°C at 150-300xg (depending on cell line)

△ CRITICAL Keep the temperature cold during all subsequent steps for the preparation of the substrate cells.

5) Remove the supernatant and resuspend each cell pellet in cold PBS.

6) Take an aliquot to count the cells and centrifuge the remainder again.

7) Remove supernatant and resuspend the cell pellet in cold freezing medium at $\sim 1 \times 10^6$ cells/mL.

8) Prepare small vials. For instance, 0.3 mL (containing approximately 330,000 cells) in 1.5 mL microtubes. Each aliquot will have enough cells for 20 gels in 2 gels/slide format. Larger aliquots can be prepared in case you plan to run more gels or slides per assay.

9) Cryopreserve at -80°C (there is no special requirement for freezing procedure; the vials can be slowly frozen using Mr Frosty® containers with isopropanol or a thick-walled polystyrene box).

10) At the day of analysis, thaw the vial and embed the cells directly in agarose (i.e. washing steps and re-suspension in new medium are not necessary).

△ CRITICAL It is essential to prepare a large batch of substrate cells (several million, depending on the number of experiments planned and the format used – see below) and they should be stored frozen in aliquots (-80°C), so that all extract samples in an experiment or a trial will be analysed on an identical substrate from the same batch. The doses of DNA-damaging agents given below are suggestions; they should be tested in dose-response experiments to confirm that the level of DNA damage is optimal for the assay (also see the Checkpoint 1).

■ PAUSE POINT Step 1 can be performed at any time before starting the experiment. The prepared substrate cells can be stored in aliquots at -80°C for at least 6 months.

Check point 1: Checking of substrate cells – day 0, • Timing ~8 h (if scoring is performed on the same day)

11) Nucleoids from Ro 19-8022 and KBrO₃ exposed cells can be checked via incubation with Fpg or hOGG1 enzymes.

12) Nucleoids from UVC-exposed cells are checked via incubation with T4 endo V enzyme.

13) Nucleoids from BPDE exposed cells can be checked with an extract of DNA repair proficient cell lines – there are currently no enzymes available that can be used to detect bulky DNA adducts.

△ CRITICAL Ideally, the lesions in the substrate should be more than enough to saturate the assay, i.e. they should result in maximum % tail DNA when incubated with the control enzyme. But on the other hand, there should not be a substantial background of SBs seen without enzyme.

△ CRITICAL Such control incubations (Steps 11-13) serve also as an *incubation reaction control* (Figure 2 – Step 4 and Figure 3), and are included in each experiment, to confirm that there really is an excess of damage for the extract to work on.

- 14) A sample of non-exposed substrate cells incubated with the corresponding enzyme – e.g. Fpg or T4endoV – will show the background level of endogenous lesions (normally insignificant in the case of T4endoV; but various cell types contain variable levels of oxidised bases) (Figure 2 – Step 4 and Figure 3).

? TROUBLESHOOTING

■ PAUSE POINT Step 11-14 can be performed at any time during the storage of substrate cells, but before starting the analysis of the extracts.

Step 2: Protein extract preparation – day 0, ● Timing ~4 h – 1 d (depending on the number of samples)

- 15) Collect biological material

(A) *Peripheral blood mononuclear cells (PBMC)*: Approximately 5-10 mL of blood is needed for preparing extract. Isolate PBMC from whole blood according to standard procedure. Approximately 1×10^6 PBMC are isolated per 1 mL of whole blood.

(B) *Cultured cells*

- (i) At least $5-10 \times 10^6$ cells are required.
- (ii) Cells grown in suspension should be in log phase of growth.
- (iii) Cells grown in monolayer in dishes or flasks should be collected by trypsinisation at sub-confluence. Centrifuge at 150-300xg for 5 min.

(C) *Solid tissues* (of animal or human origin) should be snap-frozen in liquid nitrogen as soon as possible after sampling. They can be further stored at -80°C .

- 16) Freeze and store biological material for extract preparation

(A) *Freezing cells in medium*

- (i) Isolated cells can be stored in suspension in freezing medium.
- (ii) Wash cells with PBS and count a sample. Centrifuge at 150-300xg for 5 min at 4°C .

830 (iii) Suspend cells in cold freezing medium at 5×10^6 cells/mL and prepare 1 mL aliquots in 1.5
 831 mL microtubes.

832 (iv) Freeze slowly to -80°C using a Mr Frosty® freezing container, or in a thick-walled box of
 833 expanded polystyrene.

834 *(B) Freezing cell pellets*

835 (i) Isolated cells can be snap-frozen as pellets.

836 (ii) Suspend in PBS and count a sample. Centrifuge at 150-300xg for 5 min at 4°C .

837 (iii) Suspend cells in cold PBS at 5×10^6 cells/mL and prepare 1 mL aliquots in 1.5 mL microtubes.

838 (iv) Centrifuge at $\sim 2,000 \times g$ for 5 min at 4°C .

839 (v) Carefully remove as much supernatant as possible, without disturbing the pellet.

840 (vi) Drop the tubes into liquid nitrogen. Store at -80°C . Alternatively the tubes can be placed
 841 directly at -80°C .

842 *(C) Grinding of solid tissues*

843 (i) At a later point, using a pestle and mortar pre-chilled in liquid nitrogen, grind the frozen
 844 tissue under liquid nitrogen, and divide into aliquots of about 30-50 mg using a chilled
 845 spatula.

846 **! CAUTION** Wear protective glasses and gloves when grinding tissues under liquid nitrogen.

847 (ii) Quickly weigh the frozen aliquots. Take care not to thaw them!

848 (iii) Store ground tissue aliquots at -80°C .

849 **Δ CRITICAL** Don't add any storage solution to the tissue for freezing.

850 **Δ CRITICAL** If amounts of tissue are already < 50 mg - for example, small rodent organs such as
 851 hippocampus, and tissue biopsies - do not grind, since too much of the material will be lost in
 852 the mortar. Half a hippocampus or a biopsy is typically around 5 mg, which is enough to run
 853 two assays.

854

855 **■ PAUSE POINT** Steps 15-16 can be performed at any time before starting the experiment. The
 856 harvested biological material can be either stored at -80°C for several months or used immediately,
 857 i.e. proceed to Step 17.

858

859 17) Protein extraction

860 **Δ CRITICAL** Keep samples on ice during the whole procedure of extract preparation

861 *(A) From frozen cells in freezing medium (either PBMC or cultured cells)*

862 (i) Thaw the frozen cells (either at 37°C or at room temperature), and as soon as the ice has
 863 melted, centrifuge at 150-300xg for 5 min at 4°C .

- 864 (ii) Resuspend the pellet in cold PBS and spin again.
- 865 (iii) Resuspend once more in cold PBS and centrifuge at ~2,000xg for 5 min at 4°C.
- 866 (iv) Discard supernatant by tipping out, and carefully remove the last microlitres with a pipettor
- 867 without disturbing the pellet. The pellet should be almost dry.
- 868 (v) Proceed to the next step.

869 (B) From fresh or frozen cell pellets

- 870 (i) Add 50 µl of buffer A to each pellet of 5x10⁶ cells.
- 871 (ii) Vortex vigorously, snap-freeze by dropping into liquid nitrogen and immediately thaw
- 872 again.
- 873 (iii) To each 50 µL aliquot, add 15 µL of buffer A/1% Triton® X-100.
- 874 (iv) Vortex for 5 sec and leave for 10 min on ice.
- 875 (v) Centrifuge at ~15,000xg for 5 min at 4°C to remove cell debris.
- 876 (vi) Collect the supernatant in a new microtube.
- 877 (vii) Store at -80°C or proceed to Step 4.

878 (C) From snap frozen/frozen ground tissues

- 879 (i) Thaw tissue aliquots.
- 880 (ii) Add 100 µL of buffer A per 30 mg.
- 881 (iii) Vortex vigorously, snap-freeze by dropping into liquid nitrogen and immediately thaw
- 882 again.
- 883 (iv) To each 100 µL aliquot, add 30 µL of buffer A/1% Triton® X-100.
- 884 (v) Vortex vigorously and incubate for 10 min on ice.
- 885 (vi) If necessary, larger particles of tissue can be homogenized with a microtube pestle.
- 886 (vii) Centrifuge at ~15,000xg for 5 min at 4°C to remove cell debris.
- 887 (viii) Collect the supernatant in a new microtube.
- 888 (ix) Store at -80°C or proceed to Step 4.
- 889

890 For preparation of protein extracts from *Drosophila melanogaster*, see the Box 3. A great advantage

891 of using this model is that several efficient and deficient strains are available for most of the DNA repair

892 pathways.

893 ■ **PAUSE POINT** Step 17 can be performed at any time before starting the experiment. The prepared

894 sample extracts can be either stored at -80°C for several months or used immediately, i.e. proceed to

895 the Step 18.

△ CRITICAL In case there is a large number of samples in one experiment/trial, from which the extracts are not possible to prepare in one day, it is essential to store all of them until later use. Tip: In case the protein concentration will be quantified later (step 18), store a small aliquot of the extract separately to avoid extra freeze-thaw cycles of the main extract.

? TROUBLESHOOTING

Checkpoint 2: Measure protein concentration – day 0, ● Timing ~2 h – 1 d (depending on the number of samples)

18) Quantify protein concentration

Although we advise the Lowry-based BIO-RAD DC Protein Assay Kit, using BSA as a standard, other assays might be as suitable to quantify the protein concentrations.

△ CRITICAL When selecting a protein quantification assay, choose one that: i) controls for the presence of Triton® X-100, DTT, and EDTA (detergents) in the buffer; and ii) which allows you to measure the protein concentrations at 650–750 nm to avoid interference of haemoglobin and bilirubin.

△ CRITICAL When storing the extracts for later use, it is advisable to measure protein concentration after the extract preparation, but to store it undiluted. Dilution with buffer should be performed just before the reaction.

■ PAUSE POINT Step 18 can be performed at any time before the evaluation of the DNA repair activity of the extracts.

△ CRITICAL Protein extracts can be kept on ice for immediate use or stored as smaller aliquots (to reduce freezing/thawing cycles) at -80°C for later use in the incubation reaction. Extracts can be stored at -80°C. Extracts have been stored for 1.5 months at -80°C without losing enzyme activity^{12,18}. Longer periods may also be appropriate but need to be tested.

△ CRITICAL Tissue extracts with a protein concentration of ~20-30 mg/ml can be obtained from ~50 mg of tissue, which is sufficient material for several assays (approximately 20 assays when running samples in duplicate in a 2 gels/slide format).

Step 3: Embedding substrate cells in LMP agarose & cell lysis – day 1, ● Timing ~4 h

Prepare materials

- 19) Submerge the required number of LMP agarose aliquot in boiled water to melt the agarose and then cool to 37°C (in water bath or thermoblock).
- 20) Cool centrifuge to 4°C.
- 21) Prepare working lysis solution (100 mL is needed for Coplin jar of 16 slides): to 99 mL lysis stock solution (4°C) add 1 mL Triton® X-100, mix, put into a Coplin jar, store at 4°C until use.
- 22) Put metal plate on box with ice.
- 23) Label the required number of slides on the frosted end using a pencil, not a pen.

Embedding cells in agarose and lysis

- 24) Thaw an aliquot of frozen substrate cells, both non-exposed and exposed.
- 25) As soon as the aliquot is thawed, add 1 mL of cold PBS to the 1.5mL microtube and spin at 150-300xg for 5 min at 4°C to wash cells.
- 26) Suspend pellets in cold PBS and spin again.
- 27) Remove the supernatant, disperse the pelleted cells by tapping vigorously, and add required volume of 0.7% LMP agarose (dissolved in PBS) at 37°C to reach the concentration of 2×10^5 cells/mL. For an aliquot of 3×10^5 cells, simply add 1.5 mL of agarose. Mix by pipetting gently up and down once. *Note:* Alternatively, take 45 µL of cell suspension and mix it with 105 µL of 1% low melting point agarose at 37°C (see video <https://youtu.be/T42JOvD2MnE>). This option is often applied when working with a large number of samples, so that cells can be kept on ice until use.
- 28) From each LMP-cell suspension (containing non-exposed or exposed cells), transfer two 70 µL drops to each pre-coated microscope slide (the final number of cells per gel is ~14,000).
- 29) Cover gels with 20x20mm coverslips and keep for 5-10 min at 4°C.
- 30) Remove coverslips and place slides for 1 h in lysis solution in a Coplin jar at 4°C.

▲ **CRITICAL** It is also critical to be quick, so the gels do not set before putting the coverslip.

■ **PAUSE POINT** The slides can be left in lysis solution for between 1 and 48 h but should be kept the same lysis duration for a whole set of experiments.

? TROUBLESHOOTING

Note: Alternatively, you can increase the throughput of the assay (for convenience in large epidemiological studies with many samples to be processed in a short time) using the assay with 12 gels/slide. The 12 mini-gels/slide format combined with the 12-Gel Comet Assay Unit⁸⁰ offers a perfect solution to process a large amount of samples. See Box 2 for instructions.

962

963 **Step 4: Incubation reaction – day 1, • Timing ~2 h**

964 ***Prepare materials***

965 31) Have ready a moist box in a 37°C incubator, containing suitable racks above water to ensure
966 humidity without the slides getting wet. Alternatively, put the slide moat at 37°C.

967 32) Dilute all sample extracts to the same optimized protein concentration using Buffer A/0.25%
968 Triton® X-100.

969 33) Dilute an aliquot of the 10x buffer B and/or buffer N stock in water to 1x working solution.

970 ***Detection of DNA incision activity of the extract***

971 34) Wash slides in buffer B or N (depending on the repair pathway to be studied), 3 changes, 5 min
972 each at 4°C (using Coplin jar).

973 35) Place slides on a metal plate on ice to prevent premature enzyme activity when the extract is
974 added.

975 36) Preparing sample extracts and control enzymes for incubation reaction. For a 2 gels/slide format
976 it is advised to prepare 250µL of extract mixed with incubation reaction buffer.

977 ***(A) To study BER***

978 (i) Add to the extract 4 volumes of incubation reaction buffer B.

979 (ii) Prepare a control solution: buffer A/0.25% Triton® X-100, mixed with buffer B in a ratio 1:4.

980 ***(B) To study NER***

981 (i) Add to the extract 4 volumes of buffer N.

982 **Δ CRITICAL** It is advised to add ATP to buffer N in a ratio 7:1, so the final working
983 concentration of ATP in the extract is 2.5 mM. This is not so important for freshly prepared
984 extract, but crucial for frozen extracts – since ATP degrades during long-term storage.

985 (ii) Prepare a control solution: buffer A/0.25% of Triton® X-100, mixed with buffer N+ATP in a
986 ratio 1:4.

987 **Δ CRITICAL** Keep diluted extracts and control solution on ice until use.

988 37) Add 50 µL of diluted sample extract or control solution to each gel (containing nucleoids of
989 either non-exposed or exposed cells; see Figure 3). Incubate duplicate aliquots of each sample
990 (i.e. two gels).

991 38) Cover with coverslips (22 x 22 mm for each gel or 24 x 60 mm to cover both gels).

992 39) Incubate at 37°C in a moist box in the incubator or slide moat for the required time. *Note:* The
993 incubation time is generally around 30 min but needs to be tested/optimized (see Material
994 setup for instructions).

995

996 For incubation reactions using 12 gels/slide format, see the Box 2.

997 For the use of *Drosophila melanogaster* protein extracts, see the Box 3.

998 To study DNA synthesis and ligation activity, an optional parallel step can be included (see Box 5).

999

1000 **Step 5: Comet formation – day 1, • Timing ~3 h (including washing steps)**

1001 ***Alkaline treatment & Electrophoresis***

1002 40) After the incubation of extract with substrate nucleoids, place slides immediately on ice to stop
1003 the enzyme reactions.

1004 41) Remove the coverslips and keep on ice until alkaline treatment. *Note:* Alternatively, microscope
1005 slides can be transferred directly to the electrophoresis tank.

1006 42) Incubate in cold electrophoresis solution for 40 min at 4°C.

1007 43) Electrophoresis at ~1V/cm for 30 min at 4°C. *Note:* Alternatively use the EPT that you estimated.

1008 **▲ CRITICAL** Voltage gradient should be measured across the platform carrying the slides. The time
1009 of electrophoresis should be tested before starting an experiment or a trial. Please see equipment
1010 setup.

1011 ***Neutralization & Washing***

1012 44) Neutralise gels by washing slides for 10 min in the neutralising solution (ice-cold PBS) and 10 min
1013 in ice-cold dH₂O at 4°C (use Coplin jar, or lay slides flat in a dish).

1014 *Note:* use EtOH for 12 gels/slide format (see the Box 2).

1015 45) Allow gels to dry overnight.

1016 **■ PAUSE POINT** The slides can be stored for a long time and stained and scored at any time point.

1017

1018 **Step 6: Comet visualization & Analysis – day 2, • Timing ~2 h – several days** 1019 **(depending on the number of samples)**

1020 ***Comet visualization***

1021 46) Slides can be stained with various dyes (see the section Reagents: Reagents for comet
1022 visualization).

1023 47) For staining with ethidium bromide (0.01 µg/ mL in water), or DAPI (1µg/ mL in water) - add 20-
1024 40 µL of staining solution to each gel and cover with a coverslip.

1025 48) For staining with SYBR® Gold or SYBR® Green, which give intense fluorescence, it is
1026 recommended to immerse slides in a bath of the dye at a dilution of 1:10,000 in TE buffer for 20
1027 min, followed by two 10 min washes with dH₂O. Slides are left to dry, and for viewing, 20 µL of
1028 dH₂O is added to each gel and covered with a cover slip. Alternatively, SYBR® Gold can also be
1029 added as 50 µL of the 1:10,000 dilution on top of each gel and cover with a coverslip.

1030 **! CAUTION** All dyes may be mutagenic or even carcinogenic.

1031 49) Subsequently, comets are visualized with a fluorescence microscope.

1032 **Comet analysis**

1033 50) Computer-assisted image analysis, using commercially available software, gives the most
1034 reproducible results. We recommend the use of tail intensity (TI), representing % of DNA in the
1035 tail of the comet, as the parameter to describe the comets.

1036 51) Score at least 50 comets per gel, i.e. 100 comets per sample when working in duplicates.

1037 52) To obtain the TI value per sample or control, calculate first the median TI for each gel over the
1038 scored comets (i.e. the 50 comets in each gel) and then the mean TI over the replicate gels.

1039 **△ CRITICAL** In an experiment or a trial, all comets should be scored by the same person to minimize
1040 inter-operator variation using the same software all the experiment/trial.

1041

1042 **■ PAUSE POINT** Slides can be stained and scored on the day of the experiment (day 1) or stored
1043 (un)stained.

1044

1045 **BOX 2 - 12 GELS/SLIDE FORMAT**

1046 **12 mini-gels/slide format:** 12 mini-gels per slide (in 2 rows of 6 gels each), instead of two large gels,
1047 containing the substrate cells are placed in a microscope slide (see Figure 7). Each of the gels contains
1048 about 200-250 cells, from which 100 comets will be analysed (if not using an automated image analysis
1049 system).

1050 **12-Gel Comet Assay Unit:** Each mini-gel can be incubated separately, with different extracts or control
1051 solutions, by isolating them using a silicon gasket, with the corresponding 12 holes, on top of the slide.
1052 A metal base with a guide on it is used to define the right position of the gels. The metal base, the slide
1053 containing the mini-gels, the silicon gasket and a top plate (with the corresponding 12 holes on it) are
1054 clamped together with metal clamps to isolate each mini-gel in a well and to avoid leakage between
1055 them (see Figure 7). A silicone seal/empty microscope slide on top is used to prevent evaporation of
1056 the solutions added to the wells.

1057 Specification for the use of the 12 mini-gels/slide and 12-Gel Comet Assay Unit:

1058 **PROCEDURE:**

1059 ***Step 1: Preparation of substrate cells***

1060 - After treatment of substrate cells, suspend them at 2.5×10^3 cells in 0.5 mL in cold freezing medium
1061 (This is enough for 20 slides containing 12 mini-gels each) and freeze them.

1062 ***Step 3: Embedding substrate cells in LMP agarose gel on microscope slides & cell lysis***

1063 ***Embedding cells in agarose and lysis***

1064 - After thawing and washing the cells, remove the supernatant and add 0.7% low melting point agarose
1065 to a final concentration of 0.5×10^5 cells/mL. Place an agarose-precoated slide in the ice-cold metal base
1066 and, using a multi-dispensing pipettor, place 12 drops of 5 μ L of cell suspension (containing 250 cells
1067 approximately) following the guide (2 rows of 6 gels each). Keep for 2-3 min at 4°C.

1068 *Note:* No coverslip is used; the mini-gels have domed shape.

1069 *Note:* The two rows can contain exposed or non-exposed substrate cells, or alternatively, one row of
1070 exposed and one row of non-exposed substrate cells.

1071 ***Step 4: Incubation reaction - Detection of DNA incision activity of the extract***

-Place slide/s in the 12-gel comet assay unit/s, on ice. Add 30 µL of extracts or control solutions in the corresponding wells. Cover the chamber/ with the silicone seal/ empty microscope slide and place it/them in an incubator at 37°C for the required time. The incubation time is normally around 10-30 min. After then, put the chamber/s on ice, remove the slide/s containing the gel from it/them and immediately transfer it/them to cold electrophoresis solution (on the electrophoresis tank or in Coplin jars).

Note: results from enzyme titration are different when using 2 gels/slide or 12 mini-gel/slide plus the 12-Gel Comet Assay Unit⁷⁹. Perform the titration experiments using the gel format and the equipment that is going to be used during the analysis of the samples.

Note: To stop the reaction, and even more important, to avoid cross-contamination, slides should be quickly removed from the 12-Gel Comet Assay Unit and rinsed in the electrophoresis solution before transferring slides to the electrophoresis tank (for the alkaline treatment & electrophoresis step).

Step 5: Comet formation (including washing steps)

Neutralization & Washing

- After neutralization and washing, mini-gels should be dehydrated for 15 min in 70% ethanol followed by 15 min in absolute ethanol.

? TROUBLESHOOTING

1091

1092 **BOX 3 – APPLICATION OF THE COMET-BASED *IN VITRO* DNA REPAIR ASSAY ON**
1093 ***DROSOPHILA MELANOGASTER***

1094 Specification for the use of *Drosophila melanogaster*⁸¹:

1095 **Step 2: Extract preparation and Step 4: Incubation reaction**

1096 - Place about 150 adult flies (between 7 and 12 days after eclosion/hatching), previously
1097 anesthetized, in a cold mortar on ice.

1098 *Note:* the use of a mix of males and females avoids any effect of sex.

1099 - Add 500 µL of buffer A adjusted at pH 8 and smash the flies with a pestle keeping the mortar on ice.

1100 - Split the obtained solution in 10 aliquots of 50 µL in cryotubes (discard the solid part) and place
1101 them in liquid nitrogen immediately. Store the aliquots at -80 °C until the day of the analysis.

1102 On the day of the analysis:

1103 - Place on ice an enough number of aliquots (extract obtained from one aliquot is enough for the
1104 treatment of 8 big gels, 4 slides) for thawing.

1105 - Add 12 µL of 1% Triton® X-100 per aliquot, vortex for 5 s, and place on ice for 5 min.

1106 - Centrifuge at about 15,000 x g for 10 min to separate cell debris.

1107 - Remove 50 µL supernatant and add it to 200 µL of cold Buffer B (for BER) or Buffer N (for NER).
1108 Keep the tubes on ice.

1109 *Note:* The protein concentration should be around 2.5 µg/µL.

1110 - Add 30 µL of extracts or control solutions on top of each gel containing the substrate cells, cover
1111 them with coverslips (22x22 mm), place the slides on a humidity chamber and incubate them at 24°C
1112 for 30 min.

1113

1114

1115
1116
1117
1118
1119
1120
1121
1122
1123
1124
1125
1126
1127
1128
1129
1130
1131
1132
1133
1134
1135
1136
1137
1138
1139
1140
1141
1142

BOX 4: LIST OF CELL TYPES OR CELL LINES PREVIOUSLY USED AS SUBSTRATE CELLS BY THE AUTHORS

Note: Any cell type can be used, but it is advisable to use cells in suspension to avoid trypsinization and centrifugation steps.

- **WBC** fractions can be used, but keep in mind to isolate enough cells to be able to create a large enough batch of substrate cells for the whole series of experiments. **PBMC** and **lymphocytes** have both been used successfully to assess BER activity, using Ro 19-8022 + light to induce oxidized DNA bases^{19,42,66,82}.

Cell cultures growing in suspension:

- **THP-1** cells have been exposed to KBrO₃ to induce DNA oxidation damage^{9-11,83} and UVC to assess NER (unpublished data).

- **TK6** cells have been exposed to UVC to assess NER activity, treated with MMS to assess repair of DNA alkylation lesions, and can also be treated with Ro19-8022 + light^{19,84}.

Attaching cells:

- **HeLa** cells: (i) treated with photosensitizer Ro 19-8022 + light, were used in several studies to assess BER activity; and (ii) exposed to UVC to study NER activity^{7,13,18,23,34-36,39,40,64,69,85-89}

- **A549** cells were successfully used as substrate cells using BPDE^{12,90}, and UVC (unpublished data) as DNA damage inducing agents.

- **HepG2** cells have been used as substrate cells upon UVC-exposure and treated with Ro19-8022 + light⁹¹.

- **HCT116** cells were treated by photosensitizer Ro 19-8022 + visible light to assess BER⁴⁶.

- **Caco-2** were exposed to paraquat to assess the repair of DNA oxidation lesions⁹².

1143

1144 **BOX 5: DETECTION OF DNA SYNTHESIS & LIGATION ACTIVITY OF THE EXTRACT**

1145 **Optional step in parallel to Step 4:**

1146 As explained in the introduction, assessing the DNA repair incision activity is already sufficient to study
1147 the effect of external and internal factors on an organism's DNA repair activity. For instance, a
1148 significant correlation was observed between BPDE–DNA adduct removal as studied by ³²P-post-
1149 labelling (representing the whole DNA repair process from DNA incision to ligation) and the NER
1150 incision activity as measured by the comet-based *in vitro* DNA repair assay¹². These data confirm that
1151 the DNA incision activity detected by the comet-based *in vitro* DNA repair assay is representative for
1152 an individual's DNA repair capacity.

1153 Still, if one desires to show that DNA synthesis & ligation can take place, a parallel incubation of extracts
1154 supplemented with dNTPs can be performed⁶. Note that in this case a parallel set of slides should be
1155 prepared.

1156 - First incubate this parallel set of slides with the same sample extract-reaction buffer mixture, as
1157 explained in step 35-39, for DNA incisions to occur and accumulate.

1158 - After the first incubation (of 30 min – or as optimized), rinse the slides quickly in reaction buffer (B or
1159 N, depending on the type of assay).

1160 - Add dATP, dCTP, dGTP and dTTP at 5 μM each and 2.5 mM ATP to the sample extract-reaction buffer
1161 mixtures, prepared as described in Step 36.

1162 - Put 50μL of this mixture on each gel and cover with coverslips (22 x 22 mm for each gel or 24 x 60
1163 mm to cover both gels).

1164 - Incubate at 37°C in the incubator + moist box / slide moat for another 30 min (or as optimized). So,
1165 the total incubation time will then be 60 min: 30 min incubation with extract mixtures without dNTPs
1166 + ATP, followed by 30 min incubation with extract mixture including dNTPs + ATP.

1167

1168 In case the detected incisions/SBs are back to background levels, this indicates that DNA synthesis &
1169 ligation occurred efficiently.

1170

1171 **TIMING**

1172 **Day 0:**

- 1173 Step 1: Preparation of substrate cells: 1-10, 2-6 h (depending on the DNA damaging agent)
- 1174 Checkpoint 1: Checking of substrate cells: 11-14, 8 h (if scoring is performed the same day)
- 1175 Step 2: Protein extract preparation: 15-17, 4 h - 1 d (depending on the number of samples)
- 1176 Checkpoint 2: Measure protein concentration: 18, 2 h – 1 d (depending on the number of samples)
- 1177

1178 **Day 1:**

- 1179 Step 3: Embedding substrate cells in LMP agarose gel & cell lysis: 19-30, 4 h
- 1180 Step 4: Incubation reaction: 31-39, 2 h
- 1181 Step 5: Comet formation: 40-45, 3 h
- 1182

1183 **Day 2:**

- 1184 Step 6: Comet visualization & Analysis: 46-52, 2 h - several days (depending on the number of
- 1185 samples)
- 1186
- 1187

TROUBLESHOOTING

Step	Problem	Possibly reason	Solution
Material setup: Pre-coating microscope slides; 1-4	Agarose does not hold on the slides	The presence of grease or dust on the slides	Degrease the slides by washing them with EtOH. Leave them to dry at room temperature or pass the slide through the flame of a Bunsen burner.
Step 1: Preparation of substrate cells; 2 (C)	Low levels of UVC-induced lesions (detected by the use of T4endoV)	Use of a cell line with low sensitivity to UVC (e.g., A549 cells are less sensitive than TK-6 and THP1 cells)	Increase the dose of UVC (e.g., for A549 cells 2 J/m ² was needed instead of 1 J/m ²).
Checkpoint 1: Checking of substrate cells; 11-14	Excessive level of background lesions	Use of old cells (i.e., with a high number of passages), problems in the exposure, or to high concentration of lesion-specific enzyme	Prepare another batch and/or use cells with a lower number of passages. Make sure to titrate the lesion-specific enzyme to determine the optimal conditions.
Step 3: Embedding substrate cells in LMP agarose & cell lysis; 24-30	Losing gels while removing the coverslip	Gels do not set properly due to air condensation in rooms with high temperature and/or humidity	Cool down the working room ideally to about 20°C. Embedding substrate cells in gels in an air-conditioned room is a good option.
Step 6: Comet visualization & Analysis; 50-51	Lack of incision activity of protein extracts	Presence of proteases in the sample extract	Use protease inhibitors during sample extract preparation or keep samples on ice during the whole procedure to minimise the activity of proteases in the extracts.
Box 2	Comet tails are going in all directions in the 12 mini-gels (referred to as “edge effect”)	Uneven drying of the mini-gels	Take care to dry the gels using EtOH immediately after the neutralisation. Dehydration is crucial to avoid the edge effect ⁹³ .

ANTICIPATED RESULTS:

Data analysis

Calculation of repair rates

There are 4 TI (%) values needed to calculate repair rate:

- (1) The TI (%) for the background control, which represents background DNA breaks (Microscope slide 1 – in Figure 8);
- (2) The TI (%) for the treatment control, which represents background DNA breaks plus breaks induced by the exposure to the DNA damage-inducing agent (Microscope slide 2 – in Figure 8);
- (3) The TI (%) for the specificity control, which represents background DNA breaks plus non-specific incisions induced by the sample extract (Microscope slide 3 – in Figure 8);
- (4) The TI (%) for the sample extract incubated with exposed substrate, which represents background DNA breaks, exposure-induced breaks, non-specific breaks plus specific extract incisions at lesion sites (Microscope slide 4 – in Figure 8).

To calculate the DNA repair activity, first subtract the background control value (1) from all other data (2, 3, 4), giving 'net' breaks. Subsequently, calculate repair rate by:

DNA repair incision activity =

net TI (%) for slide (4) – net TI (%) for slide (3) – net TI (%) for slide (2)

Example calculation (based on Figure 8)

The TI of the background control (1) is usually low (e.g. 2%), whereas the highest TI is usually found in the sample extract incubated with exposed substrate (4) (e.g., 23%). Reaction 2 and 3 usually yields slightly increased values when compared to incubation (1) (e.g., 5% and 3.5%). The repair rate is in this particular case: $\Delta\text{TI} = (23-2) - (5-2) - (3.5-2) = 21 - 3 - 1.5 = 16.5\%$. Results are normally reported as change in TI in a given time.

Negative values could be generated when calculating the 'net' breaks (or incisions). This can be due to the high levels of background DNA breaks in background controls. If values are close to 0 (between 0 and -2), they can be adjusted to 0 for the calculation of 'DNA repair incision activity'. However, if values are lower, the experiment should be repeated just in case technical problems have occurred. If results are the same, excessive level of background lesions can be present (see 'Troubleshooting table, Checkpoint 1: Checking of substrate cells; 11-14'). If the calculated DNA repair incision activity of an

extract is negative, but all the internal experimental controls give the expected results and the optimization of protein extract and time of incubation is correct, it can be assumed that the extract does not have DNA repair incision activity or that it is below the limit of detection of this assay.

Example data

Figure 9 depicts example data of a successful (Figure 9A) and a suboptimal (Figure 9B) assay that assessed the BER activity of piglet hippocampus samples, using substrate cells (HeLa) exposed to Ro 19-8022 + light or non-exposed substrate cells (PBS + light)³⁰. The problem with the assay results shown in Figure 9B involves the high values observed for the non-exposed substrate cells. Since the incubation of exposed substrate cells with Fpg (incubation reaction control) showed expected results (% tail DNA = ~60-70%), it is suggested that there was not a technical error but possibly an error in the handling of the substrate cells or preparation of those slides.

Since KBrO₃ induces the same type of lesions as Ro 19-8022 + light, similar data can be expected when using KBrO₃-exposed substrate cells. The use of KBrO₃ in the comet-based *in vitro* DNA repair assay is increasing, and example data can be consulted in various reports showing DNA repair activity in cell cultures⁹, animal tissues¹⁰ and human PBMC¹¹.

Example data for NER assays are shown in **Figure 10** for the use of UVC-exposed substrate cells and **Figure 11** for the use of BPDE-exposed substrate cells. Figure 10A shows results from PBMC extracts from 7 volunteers and all assay controls indicate that the assay performed well. In contrast, Figure 10B shows results from saliva cell extracts from the same 7 volunteers. Here, the results of the incubation with T4endoV indicate that the incubation reaction performed well, but the background and treatment control are higher than expected, which could be due to the handling of the cells. Regardless, the specificity controls indicate that most extracts suffer from the presence of non-specific nuclease activity – which in the case of saliva could be due to the presence of bacterial enzymes. This assay would need further optimization to reduce this non-specific incision activity.

In **Figure 11** anticipated results from PBMC extracts incubated with BPDE-exposed substrate cells are shown, which illustrate the importance of checking the substrate cells (as described in checkpoint 1). When A549 cells were treated with 0.5% DMSO for 30 min at 37°C to create non-exposed substrate cells, increased background levels were observed – probably due to the toxicity of DMSO (Figure 11B). While treatment with 0.12% DMSO for 30 min at 4°C, in parallel to exposure to 1µM BPDE for 30min at 4°C, resulted in suitable substrate cells with low background levels (Figure 11A). Previously, PBMC samples from these 8 volunteers were tested, and the one showing the highest DNA repair incision activity was selected to serve as incubation reaction control (IC) for that series of experiments.

66

67 ***How to express results?***

68 The incubation of the sample extract with substrate DNA should be performed for a set length of time
69 within the same series of experiments (optimised as described in “Material setup”). So, the results can
70 be expressed as a rate of accumulation of DNA SBs which are optimally expressed as % of DNA in tail.
71 Comet assay results can be further converted to an actual DNA break frequency, using a calibration
72 curve based on irradiated cells, so that the results can be expressed in terms of breaks per 10^9 Da^{94,95}.
73 Another way of expressing results is to calculate DNA SBs relative to protein concentration, e.g.
74 amount/number of breaks per mg/mL protein. Since the same volume of the sample extract is added
75 to each gel containing the approximately the same amount of substrate DNA, the incision activity can
76 also be expressed as an amount/number of breaks per mg protein. It is, however, important to keep
77 in mind that activity, as assayed with this method, is not linearly proportional to protein concentration
78 but increases less than 2-fold for a doubling of concentration^{18,67,96}.

79 The latter point makes normalisation of the results difficult, but it is possible to allow correction for
80 inter-assay variations. Including the “incubation reaction controls” in each assay allows for
81 normalisation of the data. An alternative way to normalise the data is to include a sample extract from
82 a pooled/reference sample in each assay that has a protein concentration similar to the tested samples
83 and of which the detected incision activity is situated in the linear part of the titration curve. Based on
84 the data for this pooled sample, or the incubation reaction control, it is possible to correct the values
85 of the other samples as follows:

- 86 • First calculate the DNA repair incision activity, as indicated above.
- 87 • List the DNA repair incision activity for the pooled/reference sample (or the incubation reaction
88 control) as detected in all the assays within one set of experiments.
- 89 • Calculate the median value, M
- 90 • The detected DNA repair incision activity for the pooled/reference sample (or the incubation reaction
91 control) in a particular experiment X, defined as Q, will then serve in the correction factor - being M/Q.
- 92 • Multiply the values of DNA repair incision activity for each sample in experiment X by M/Q to
93 normalize all values.
- 94 • Example data and calculation can be consulted in the supplementary file 1.

95

96 The variations to be expected in the levels of detected DNA repair activity, and the effect of intrinsic
97 and external factors on DNA repair activity levels, have been described in a recent comprehensive
98 review³⁸.

99

100 **Reporting Summary**

101 Further information on research design and ethical approvals is available in the Nature Research
102 Reporting Summary linked to this article.

103

104 **Data availability**

105 The authors declare that the majority of the data shown here as examples or anticipated results are
106 available in the original papers. Other supporting data are available upon reasonable request to the
107 corresponding author. For instance, figures 10 and 11 are based on unpublished data, generated by
108 WG5 within the hCOMET-COST CA15132 action].

109

110

FIGURES:

Figure 1. Principle of the sample extract incubation reaction (for BER acting against DNA oxidation damage as example).

Abbreviations: 8-oxoG – 8-oxoguanine, AP site – apurinic/aprimidinic site.

Figure 2. Stepwise overview of the technique (using 2 gels/slide as example).

Abbreviations: dNTPs – deoxyribonucleotide triphosphates, WBC – white blood cells.

Figure 3. Overview of potential assay setup for the comet based *in vitro* DNA repair assay. This overview describes assay setup to assess BER and NER in 3 samples, using Ro 19-8022 + Light (Ro) or UVC (UV) exposed and corresponding non-exposed (noRo or nUV) substrate cells as example.

Figure 4. Equipment setup for exposure of substrate cells. Setup of the 500 Watt lamp at a distance of 33cm above the cells that are placed on an ice box, to perform the exposure to the photosensitiser Ro 19-8022 (left). Construction of a cardboard box (black) with cotton gauze to reduce the intensity during the UVC exposure.

Figure 5. Example data illustrating how to select the optimal incubation time. The red dotted line indicates that for experiments where a mouse hippocampus tissue extract of 3mg/mL is incubated with substrate cells (HeLa) exposed to 1 μ M Ro 19-8022 + light an incubation time of 25min will be the optimal timing, allowing the detection of the DNA repair incision activity in the linear part of the activity-with-time curve.

Figure 6. Principle of optimization of the extract's protein concentration. (A) A too low protein concentrations will yield low background levels, but also low DNA repair specific incisions; (B) the "optimal" protein concentration results in low background levels and high levels of DNA repair incisions; (C) at high concentrations, non-specific nucleases may increase the background.

Figure 7. The 12 gels/slide format. (A) 12-Gel Comet Assay Unit (NorGenoTech) to incubate the gels separately. From bottom: metal position guide, microscope slide, silicone gasket and top plate

142 tightened by metal clamps. Down right: twelve mini-gels set on a microscope slide; (B) Schematic
143 overview of 12 gels/slide format. Gels containing non-exposed and exposed substrate cells on different
144 slides. Incubation positions for buffer (B), sample extracts (1-10) and incubation reaction control (IC).

146 **Figure 8. Overview of slides containing internal experimental controls and sample extract used to**
147 **calculate the final DNA repair incision activity.** 1, 2, and 3 represent internal experimental controls
148 and 4 represents sample extract. The bars on the graph shows example data with SD for the internal
149 experimental controls (1, 2 and 3) and a piglet (*Sus scrofa domestica*) hippocampus protein extract
150 (4)(3mg/mL). Data are shown as mean TI of two independent incubations (i.e., on 2 separate gels)
151 within one experiment.

153 **Figure 9. Examples of results obtained with a BER assay.** Results were produced using substrate cells
154 (HeLa) exposed to Ro 19-8022 + light (grey bars) or non-exposed substrate cells (white bars) incubated
155 with various protein extracts from piglet hippocampus samples (n=5, *Sus scrofa domestica*)³⁰. (A)
156 Example of a successful experiment and (B) example of a suboptimal experiment – showing too high
157 background values. Data are shown as mean % tail DNA of two independent incubations with SD within
158 one experiment.

160 **Figure 10. Examples of results obtained with a NER assay.** Results were produced using substrate cells
161 (TK6) exposed to 1 J/m² UV (grey bars) or non-exposed substrate cells (white bars) incubated with
162 various human protein extracts either (A) from PBMC, or (B) from saliva cells. Data are shown as mean
163 % tail DNA of two independent incubations with SD within one experiment.

165 **Figure 11. Examples of results obtained with a NER assay.** Results were produced using substrate cells
166 (A549) exposed to 1μM BPDE (grey bars) or non-exposed substrate cells (white bars; exposed to vehicle
167 control DMSO) incubated with various protein extracts from human PBMC. (A) Non-exposed substrate
168 cells were treated with 0.12% DMSO for 30 min at 4°C, while (B) treatment with 0.5% DMSO for 30
169 min at 37°C increased background levels. Data are shown as mean % tail DNA of two independent
170 incubations with SD within one experiment. IC = Incubation Reaction Control.

173 **TABLES:**

174 **Table 1.** Overview of the various modifications of the comet assay.

	Standard comet assay	Enzyme-modified comet assay	Cellular repair assay	Comet-based <i>in vitro</i> DNA repair assay
Starting material to embed in the gel	Cells from samples	Cells from samples	Cells from samples	Substrate cells containing specific DNA lesions
Enzyme treatment during incubation step	No incubation step needed	DNA endonucleases isolated from bacteria or human cells	No incubation step needed, but cells are grown and harvested at various time points	DNA repair enzymes present in protein extracts from samples will incise at lesions in substrate DNA
Endpoint	Single and double SBs	Specific lesions; e.g. oxidized, alkylated, methylated, or dimerised bases	Induction and removal of DNA lesions over time.	DNA repair incision activity
Data interpretation	More migration of DNA to the tail indicates higher levels of DNA damage	More migration of DNA to the tail compared to controls (not incubated with enzyme) indicates higher levels of specific DNA lesions	More migration of DNA to the tail indicates higher levels of DNA damage. Over time, less DNA in the tail will indicate repair of SBs.	More migration of DNA to the tail compared to controls (not incubated with protein extract) indicates higher DNA repair activity

175

176 **Table 2.** Recommended protein concentrations for various cell types and tissues.

Cell type/Tissue	Protein concentration (mg/mL)	References
Lymphocytes	2	12
Fibroblasts	0.1	12
Human colon biopsies	3	19,46,69
Mouse brain & liver	3 - 5	18
Pig brain	3 - 5	30
Mouse & pig colon	0.3 - 0.5	69
Mouse lung	0.25-2.5	69,97

177

178 **Table 3.** Incubation reaction controls that can be used for checking the amount of DNA damage in
179 substrate cells and to normalize between batches of experiments.

Type of exposure	Type of DNA lesion	Potential Control enzyme
Ro 19-8022	8-oxoG	Fpg, hOGG1
KBrO ₃	8-oxoG	Fpg, hOGG1
UVC light	Thymidine dimers	T4endoV
BPDE	BPDE-G	Standard extract of DNA repair proficient cell line

180

181 **SUPPLEMENTARY INFORMATION:**

182 Supplementary Results 1 – Examples of data normalisation

183

184

AUTHOR CONTRIBUTIONS STATEMENTS:

S.V., A.A., R.G. and S.L. designed figures; S.L. provided anticipated results and managed the manuscript preparation; S.V., A.A. and S.L. drafted the paper, and A.C., M.D., P.M. A.O., I.G., and P.V. contributed to and revised the manuscript. All authors read and approved the final manuscript.

ORCID FOR CORRESPONDING AUTHORS:

ORCID: [0000-0003-3288-7331](https://orcid.org/0000-0003-3288-7331)

ACKNOWLEDGMENTS:

We thank the hCOMET project (COST Action, CA 15132) for support. A.A. thanks the Ministry of Economy, Industry and Competitiveness ('Ramón y Cajal' programme, RYC2013-14370) of the Spanish Government for personal support. P.V. acknowledges the support from National Science Foundation (19-10543S).

COMPETING INTERESTS:

The authors declare no competing interests.

202 REFERENCES:

- 203 1 Ostling, O. & Johanson, K. J. Microelectrophoretic study of radiation-induced DNA
204 damages in individual mammalian cells. *Biochem Biophys Res Commun* **123**, 291-298,
205 doi:10.1016/0006-291x(84)90411-x (1984).
- 206 2 Cook, P. R., Brazell, I. A. & Jost, E. Characterization of nuclear structures containing
207 superhelical DNA. *J Cell Sci* **22**, 303-324 (1976).
- 208 3 Collins, A. R. & Azqueta, A. Single cell gel electrophoresis combined with lesion-
209 specific enzymes to measure oxidative damage to DNA. *Laboratory Methods in Cell*
210 *Biology* **112**, 69-92, doi:10.1016/b978-0-12-405914-6.00004-4 (2012).
- 211 4 Stefanini, M. *et al.* Novel Chinese hamster ultraviolet-sensitive mutants for excision
212 repair form complementation groups 9 and 10. *Cancer Res* **51**, 3965-3971 (1991).
- 213 5 Lorenzo, Y. *et al.* The carotenoid beta-cryptoxanthin stimulates the repair of DNA
214 oxidation damage in addition to acting as an antioxidant in human cells.
215 *Carcinogenesis* **30**, 308-314, doi:10.1093/carcin/bgn270 (2009).
- 216 6 Collins, A. R., Fleming, I. M. & Gedik, C. M. In vitro repair of oxidative and ultraviolet-
217 induced DNA damage in supercoiled nucleoid DNA by human cell extract. *Biochim*
218 *Biophys Acta* **1219**, 724-727, doi:10.1016/0167-4781(94)90236-4 (1994).
- 219 7 Collins, A. R. *et al.* Inter-individual differences in repair of DNA base oxidation,
220 measured in vitro with the comet assay. *Mutagenesis* **16**, 297-301,
221 doi:10.1093/mutage/16.4.297 (2001).
- 222 8 Moller, P. *et al.* Searching for assay controls for the Fpg- and hOGG1-modified comet
223 assay. *Mutagenesis* **33**, 9-19, doi:10.1093/mutage/gex015 (2018).
- 224 9 Borghini, A., Roursgaard, M., Andreassi, M. G., Kermanizadeh, A. & Moller, P. Repair
225 activity of oxidatively damaged DNA and telomere length in human lung epithelial
226 cells after exposure to multi-walled carbon nanotubes. *Mutagenesis* **32**, 173-180,
227 doi:10.1093/mutage/gew036 (2017).
- 228 10 Jensen, D. M. *et al.* Telomere length and genotoxicity in the lung of rats following
229 intragastric exposure to food-grade titanium dioxide and vegetable carbon particles.
230 *Mutagenesis* **34**, 203-214, doi:10.1093/mutage/gez003 (2019).
- 231 11 Lohr, M. *et al.* Association between age and repair of oxidatively damaged DNA in
232 human peripheral blood mononuclear cells. *Mutagenesis* **30**, 695-700,
233 doi:10.1093/mutage/gev031 (2015).
- 234 12 Langie, S. A. *et al.* Development and validation of a modified comet assay to
235 phenotypically assess nucleotide excision repair. *Mutagenesis* **21**, 153-158,
236 doi:10.1093/mutage/gel013 (2006).
- 237 13 Gaivao, I., Piasek, A., Brevik, A., Shaposhnikov, S. & Collins, A. R. Comet assay-based
238 methods for measuring DNA repair in vitro; estimates of inter- and intra-individual
239 variation. *Cell Biol Toxicol* **25**, 45-52, doi:10.1007/s10565-007-9047-5 (2009).
- 240 14 Herrera, M. *et al.* Differences in repair of DNA cross-links between lymphocytes and
241 epithelial tumor cells from colon cancer patients measured in vitro with the comet
242 assay. *Clin Cancer Res* **15**, 5466-5472, doi:10.1158/1078-0432.CCR-08-3268 (2009).
- 243 15 van Dyk, E., Steenkamp, A., Koekemoer, G. & Pretorius, P. J. Hereditary tyrosinemia
244 type 1 metabolites impair DNA excision repair pathways. *Biochem Biophys Res*
245 *Commun* **401**, 32-36, doi:10.1016/j.bbrc.2010.09.002 (2010).

- 246 16 Langie, S. A. *et al.* The effect of oxidative stress on nucleotide-excision repair in colon
247 tissue of newborn piglets. *Mutat Res* **695**, 75-80,
248 doi:10.1016/j.mrgentox.2009.12.005 (2010).
- 249 17 Mikkelsen, L. *et al.* Aging and defense against generation of 8-oxo-7,8-dihydro-2'-
250 deoxyguanosine in DNA. *Free Radic Biol Med* **47**, 608-615,
251 doi:10.1016/j.freeradbiomed.2009.05.030 (2009).
- 252 18 Langie, S. A. *et al.* Measuring DNA repair incision activity of mouse tissue extracts
253 towards singlet oxygen-induced DNA damage: a comet-based in vitro repair assay.
254 *Mutagenesis* **26**, 461-471, doi:10.1093/mutage/ger005 (2011).
- 255 19 Slyskova, J. *et al.* Functional, genetic, and epigenetic aspects of base and nucleotide
256 excision repair in colorectal carcinomas. *Clin Cancer Res* **18**, 5878-5887,
257 doi:10.1158/1078-0432.CCR-12-1380 (2012).
- 258 20 Azqueta, A., Slyskova, J., Langie, S. A., O'Neill Gaivao, I. & Collins, A. Comet assay to
259 measure DNA repair: approach and applications. *Front Genet* **5**, 288,
260 doi:10.3389/fgene.2014.00288 (2014).
- 261 21 Cho, E. *et al.* Aop: 296 - Oxidative DNA damage leading to chromosomal aberrations
262 and mutations, <<https://aopwiki.org/aops/296>> (2020).
- 263 22 Silva, J. P., Gomes, A. C. & Coutinho, O. P. Oxidative DNA damage protection and
264 repair by polyphenolic compounds in PC12 cells. *Eur J Pharmacol* **601**, 50-60,
265 doi:10.1016/j.ejphar.2008.10.046 (2008).
- 266 23 Ramos, A. A., Azqueta, A., Pereira-Wilson, C. & Collins, A. R. Polyphenolic compounds
267 from Salvia species protect cellular DNA from oxidation and stimulate DNA repair in
268 cultured human cells. *J Agric Food Chem* **58**, 7465-7471, doi:10.1021/jf100082p
269 (2010).
- 270 24 Ramos, A. A., Pereira-Wilson, C. & Collins, A. R. Protective effects of ursolic acid and
271 luteolin against oxidative DNA damage include enhancement of DNA repair in Caco-2
272 cells. *Mutat Res* **692**, 6-11, doi:10.1016/j.mrfmmm.2010.07.004 (2010).
- 273 25 Azqueta, A., Costa, S., Lorenzo, Y., Bastani, N. E. & Collins, A. R. Vitamin C in cultured
274 human (HeLa) cells: lack of effect on DNA protection and repair. *Nutrients* **5**, 1200-
275 1217, doi:10.3390/nu5041200 (2013).
- 276 26 Silva, J. P., Gomes, A. C., Proenca, F. & Coutinho, O. P. Novel nitrogen compounds
277 enhance protection and repair of oxidative DNA damage in a neuronal cell model:
278 comparison with quercetin. *Chem Biol Interact* **181**, 328-337,
279 doi:10.1016/j.cbi.2009.07.024 (2009).
- 280 27 Sliwinski, T. *et al.* STI571 reduces NER activity in BCR/ABL-expressing cells. *Mutat Res*
281 **654**, 162-167, doi:10.1016/j.mrgentox.2008.06.002 (2008).
- 282 28 Folkmann, J. K. *et al.* Oxidatively damaged DNA in rats exposed by oral gavage to C60
283 fullerenes and single-walled carbon nanotubes. *Environ Health Perspect* **117**, 703-
284 708, doi:10.1289/ehp.11922 (2009).
- 285 29 Langie, S. A. *et al.* Maternal folate depletion and high-fat feeding from weaning
286 affects DNA methylation and DNA repair in brain of adult offspring. *FASEB J* **27**, 3323-
287 3334, doi:10.1096/fj.12-224121 (2013).
- 288 30 Langie, S. A. *et al.* Redox and epigenetic regulation of the APE1 gene in the
289 hippocampus of piglets: The effect of early life exposures. *DNA Repair (Amst)* **18**, 52-
290 62, doi:10.1016/j.dnarep.2014.03.011 (2014).
- 291 31 Langie, S. A. *et al.* The Ageing Brain: Effects on DNA Repair and DNA Methylation in
292 Mice. *Genes (Basel)* **8**, doi:10.3390/genes8020075 (2017).

293 32 Setayesh, T. *et al.* Impact of Weight Loss Strategies on Obesity-Induced DNA Damage.
 294 *Mol Nutr Food Res* **63**, e1900045, doi:10.1002/mnfr.201900045 (2019).

295 33 Gaivao, I. & Sierra, L. M. Drosophila comet assay: insights, uses, and future
 296 perspectives. *Front Genet* **5**, 304, doi:10.3389/fgene.2014.00304 (2014).

297 34 Dusinska, M., Dzupinkova, Z., Wsolova, L., Harrington, V. & Collins, A. R. Possible
 298 involvement of XPA in repair of oxidative DNA damage deduced from analysis of
 299 damage, repair and genotype in a human population study. *Mutagenesis* **21**, 205-
 300 211, doi:10.1093/mutage/gel016 (2006).

301 35 Slyskova, J. *et al.* Relationship between the capacity to repair 8-oxoguanine,
 302 biomarkers of genotoxicity and individual susceptibility in styrene-exposed workers.
 303 *Mutat Res* **634**, 101-111, doi:10.1016/j.mrgentox.2007.06.012 (2007).

304 36 Dusinska, M. *et al.* Are glutathione S transferases involved in DNA damage signalling?
 305 Interactions with DNA damage and repair revealed from molecular epidemiology
 306 studies. *Mutat Res* **736**, 130-137, doi:10.1016/j.mrfmmm.2012.03.003 (2012).

307 37 Staruchova, M. *et al.* Occupational exposure to mineral fibres. Biomarkers of
 308 oxidative damage and antioxidant defence and associations with DNA damage and
 309 repair. *Mutagenesis* **23**, 249-260, doi:10.1093/mutage/gen004 (2008).

310 38 Azqueta, A. *et al.* DNA repair as a human biomonitoring tool: Comet assay
 311 approaches. *Mutat Res* **781**, 71-87, doi:10.1016/j.mrrev.2019.03.002 (2019).

312 39 Dusinska, M. *et al.* Genotoxic effects of asbestos in humans. *Mutat Res* **553**, 91-102,
 313 doi:10.1016/j.mrfmmm.2004.06.027 (2004).

314 40 Dusinska, M. *et al.* Does occupational exposure to mineral fibres cause DNA or
 315 chromosome damage? *Mutat Res* **553**, 103-110, doi:10.1016/j.mrfmmm.2004.06.029
 316 (2004).

317 41 Vodicka, P. *et al.* Cytogenetic markers, DNA single-strand breaks, urinary metabolites,
 318 and DNA repair rates in styrene-exposed lamination workers. *Environ Health Perspect*
 319 **112**, 867-871, doi:10.1289/ehp.6849 (2004).

320 42 Stoyanova, E. *et al.* Base excision repair capacity in chronic renal failure patients
 321 undergoing hemodialysis treatment. *Cell Biochem Funct* **32**, 177-182,
 322 doi:10.1002/cbf.2989 (2014).

323 43 Fikrova, P. *et al.* DNA crosslinks, DNA damage and repair in peripheral blood
 324 lymphocytes of non-small cell lung cancer patients treated with platinum derivatives.
 325 *Oncol Rep* **31**, 391-396, doi:10.3892/or.2013.2805 (2014).

326 44 Slyskova, J. *et al.* Differences in nucleotide excision repair capacity between newly
 327 diagnosed colorectal cancer patients and healthy controls. *Mutagenesis* **27**, 225-232,
 328 doi:10.1093/mutage/ger088 (2012).

329 45 Slyskova, J. *et al.* Post-treatment recovery of suboptimal DNA repair capacity and
 330 gene expression levels in colorectal cancer patients. *Mol Carcinog* **54**, 769-778,
 331 doi:10.1002/mc.22141 (2015).

332 46 Vodenkova, S. *et al.* Base excision repair capacity as a determinant of prognosis and
 333 therapy response in colon cancer patients. *DNA Repair (Amst)* **72**, 77-85,
 334 doi:10.1016/j.dnarep.2018.09.006 (2018).

335 47 Slyskova, J., Langie, S. A., Collins, A. R. & Vodicka, P. Functional evaluation of DNA
 336 repair in human biopsies and their relation to other cellular biomarkers. *Front Genet*
 337 **5**, 116, doi:10.3389/fgene.2014.00116 (2014).

- 48 Singh, N. P., McCoy, M. T., Tice, R. R. & Schneider, E. L. A simple technique for quantitation of low levels of DNA damage in individual cells. *Exp Cell Res* **175**, 184-191, doi:10.1016/0014-4827(88)90265-0 (1988).
- 49 Collins, A. R., Ord, M. J. & Johnson, R. T. Correlations of DNA damage and repair with nuclear and chromosomal damage in HeLa cells caused by methylnitrosamides. *Cancer Res* **41**, 5176-5187 (1981).
- 50 Crebelli, R. *et al.* Biomonitoring of primary aluminium industry workers: detection of micronuclei and repairable DNA lesions by alkaline SCGE. *Mutat Res* **516**, 63-70, doi:10.1016/s1383-5718(02)00028-1 (2002).
- 51 Vande Loock, K., Decordier, I., Ciardelli, R., Haumont, D. & Kirsch-Volders, M. An aphidicolin-block nucleotide excision repair assay measuring DNA incision and repair capacity. *Mutagenesis* **25**, 25-32, doi:10.1093/mutage/geb039 (2010).
- 52 Lambert, B., Ringborg, U. & Skoog, L. Age-related decrease of ultraviolet light-induced DNA repair synthesis in human peripheral leukocytes. *Cancer Res* **39**, 2792-2795 (1979).
- 53 Athas, W. F., Hedayati, M. A., Matanoski, G. M., Farmer, E. R. & Grossman, L. Development and field-test validation of an assay for DNA repair in circulating human lymphocytes. *Cancer Res* **51**, 5786-5793 (1991).
- 54 Redaelli, A., Magrassi, R., Bonassi, S., Abbondandolo, A. & Frosina, G. AP endonuclease activity in humans: development of a simple assay and analysis of ten normal individuals. *Teratog Carcinog Mutagen* **18**, 17-26 (1998).
- 55 Elliott, R. M., Astley, S. B., Southon, S. & Archer, D. B. Measurement of cellular repair activities for oxidative DNA damage. *Free Radic Biol Med* **28**, 1438-1446, doi:10.1016/s0891-5849(00)00271-9 (2000).
- 56 Roldan-Arjona, T. *et al.* Molecular cloning and functional expression of a human cDNA encoding the antimutator enzyme 8-hydroxyguanine-DNA glycosylase. *Proc Natl Acad Sci U S A* **94**, 8016-8020, doi:10.1073/pnas.94.15.8016 (1997).
- 57 Sauvaigo, S. *et al.* An oligonucleotide microarray for the monitoring of repair enzyme activity toward different DNA base damage. *Anal Biochem* **333**, 182-192, doi:10.1016/j.ab.2004.06.046 (2004).
- 58 Paz-Elizur, T. *et al.* DNA repair activity for oxidative damage and risk of lung cancer. *J Natl Cancer Inst* **95**, 1312-1319, doi:10.1093/jnci/djg033 (2003).
- 59 Paz-Elizur, T. *et al.* Reduced repair of the oxidative 8-oxoguanine DNA damage and risk of head and neck cancer. *Cancer Res* **66**, 11683-11689, doi:10.1158/0008-5472.CAN-06-2294 (2006).
- 60 Paz-Elizur, T. *et al.* Development of an enzymatic DNA repair assay for molecular epidemiology studies: distribution of OGG activity in healthy individuals. *DNA Repair (Amst)* **6**, 45-60, doi:10.1016/j.dnarep.2006.08.003 (2007).
- 61 Leitner-Dagan, Y. *et al.* N-methylpurine DNA glycosylase and OGG1 DNA repair activities: opposite associations with lung cancer risk. *J Natl Cancer Inst* **104**, 1765-1769, doi:10.1093/jnci/djs445 (2012).
- 62 Leitner-Dagan, Y. *et al.* Enzymatic MPG DNA repair assays for two different oxidative DNA lesions reveal associations with increased lung cancer risk. *Carcinogenesis* **35**, 2763-2770, doi:10.1093/carcin/bgu214 (2014).
- 63 Azqueta, A. & Collins, A. R. The essential comet assay: a comprehensive guide to measuring DNA damage and repair. *Arch Toxicol* **87**, 949-968, doi:10.1007/s00204-013-1070-0 (2013).

385 64 Collins, A. R., Harrington, V., Drew, J. & Melvin, R. Nutritional modulation of DNA
386 repair in a human intervention study. *Carcinogenesis* **24**, 511-515,
387 doi:10.1093/carcin/24.3.511 (2003).

388 65 Riso, P. *et al.* DNA damage and repair activity after broccoli intake in young healthy
389 smokers. *Mutagenesis* **25**, 595-602, doi:10.1093/mutage/geq045 (2010).

390 66 Danielsen, P. H. *et al.* Oxidatively damaged DNA and its repair after experimental
391 exposure to wood smoke in healthy humans. *Mutat Res* **642**, 37-42,
392 doi:10.1016/j.mrfmmm.2008.04.001 (2008).

393 67 Abstracts of the 12th International Comet Assay Workshop held at the University of
394 Navarra, Pamplona, Spain, 29–31 August 2017 (<https://icaw.vito.be/>). *Mutagenesis*
395 **32**, e1-e28, doi:10.1093/mutage/gex037 (2018).

396 68 Guarnieri, S. *et al.* DNA repair phenotype and dietary antioxidant supplementation.
397 *Br J Nutr* **99**, 1018-1024, doi:10.1017/S0007114507842796 (2008).

398 69 Gorniak, J. P. *et al.* Tissue differences in BER-related incision activity and non-specific
399 nuclease activity as measured by the comet assay. *Mutagenesis* **28**, 673-681,
400 doi:10.1093/mutage/get047 (2013).

401 70 Jensen, A. *et al.* Influence of the OGG1 Ser326Cys polymorphism on oxidatively
402 damaged DNA and repair activity. *Free Radic Biol Med* **52**, 118-125,
403 doi:10.1016/j.freeradbiomed.2011.09.038 (2012).

404 71 Boer, K., Isenmann, S. & Deufel, T. Strong interference of hemoglobin concentration
405 on CSF total protein measurement using the trichloroacetic acid precipitation
406 method. *Clin Chem Lab Med* **45**, 112-113, doi:10.1515/CCLM.2007.013 (2007).

407 72 Roman, Y., Bomsel-Demontoy, M. C., Levrier, J., Chaste-Duvernoy, D. & Jalme, M. S.
408 Effect of hemolysis on plasma protein levels and plasma electrophoresis in birds. *J*
409 *Wildl Dis* **45**, 73-80, doi:10.7589/0090-3558-45.1.73 (2009).

410 73 Brodersen, R. Bilirubin. Solubility and interaction with albumin and phospholipid. *J*
411 *Biol Chem* **254**, 2364-2369 (1979).

412 74 Kjellin, K. G. Bilirubin compounds in the CSF. *J Neurol Sci* **13**, 161-173,
413 doi:10.1016/0022-510x(71)90145-6 (1971).

414 75 Collins, A. R. & Azqueta, A. DNA repair as a biomarker in human biomonitoring
415 studies; further applications of the comet assay. *Mutat Res* **736**, 122-129,
416 doi:10.1016/j.mrfmmm.2011.03.005 (2012).

417 76 Forchhammer, L. *et al.* Variation in assessment of oxidatively damaged DNA in
418 mononuclear blood cells by the comet assay with visual scoring. *Mutagenesis* **23**,
419 223-231, doi:10.1093/mutage/gen006 (2008).

420 77 Azqueta, A. *et al.* The influence of scoring method on variability in results obtained
421 with the comet assay. *Mutagenesis* **26**, 393-399, doi:10.1093/mutage/geq105 (2011).

422 78 Brunborg, G., Rolstadaas, L. & Gutzkow, K. B. Electrophoresis in the Comet Assay in
423 *Electrophoresis - Life Sciences Practical Applications*. (eds O-M. Boldura & C. Baltă)
424 (IntechOpen, 2018).

425 79 Muruzabal, D., Langie, S. A. S., Pourrut, B. & Azqueta, A. The enzyme-modified comet
426 assay: Enzyme incubation step in 2 vs 12-gels/slide systems. *Mutat Res* **845**, 402981,
427 doi:10.1016/j.mrgentox.2018.11.005 (2019).

428 80 Shaposhnikov, S. *et al.* Twelve-gel slide format optimised for comet assay and
429 fluorescent in situ hybridisation. *Toxicol Lett* **195**, 31-34,
430 doi:10.1016/j.toxlet.2010.02.017 (2010).

- 81 Sierra, L. M. & Gaivão, I. Chapter 23: Use of the Comet assay to study DNA repair in *Drosophila melanogaster* in *Genotoxicity and DNA repair: A practical approach*. (Humana Press, 2014).
- 82 Brauner, E. V. *et al.* Exposure to ultrafine particles from ambient air and oxidative stress-induced DNA damage. *Environ Health Perspect* **115**, 1177-1182, doi:10.1289/ehp.9984 (2007).
- 83 Moller, P. *et al.* Measurement of oxidative damage to DNA in nanomaterial exposed cells and animals. *Environ Mol Mutagen* **56**, 97-110, doi:10.1002/em.21899 (2015).
- 84 Hasplova, K. *et al.* DNA alkylation lesions and their repair in human cells: modification of the comet assay with 3-methyladenine DNA glycosylase (AlkD). *Toxicol Lett* **208**, 76-81, doi:10.1016/j.toxlet.2011.10.005 (2012).
- 85 Dusinska, M. *et al.* Testing strategies for the safety of nanoparticles used in medical applications. *Nanomedicine (Lond)* **4**, 605-607, doi:10.2217/nnm.09.47 (2009).
- 86 Choi, S. W., Yeung, V. T., Collins, A. R. & Benzie, I. F. Redox-linked effects of green tea on DNA damage and repair, and influence of microsatellite polymorphism in HMOX-1: results of a human intervention trial. *Mutagenesis* **30**, 129-137, doi:10.1093/mutage/geu022 (2015).
- 87 Brevik, A. *et al.* Supplementation of a western diet with golden kiwifruits (*Actinidia chinensis* var. 'Hort 16A') effects on biomarkers of oxidation damage and antioxidant protection. *Nutr J* **10**, 54, doi:10.1186/1475-2891-10-54 (2011).
- 88 Hanova, M. *et al.* Modulation of DNA repair capacity and mRNA expression levels of XRCC1, hOGG1 and XPC genes in styrene-exposed workers. *Toxicol Appl Pharmacol* **248**, 194-200, doi:10.1016/j.taap.2010.07.027 (2010).
- 89 Humphreys, V. *et al.* Age-related increases in DNA repair and antioxidant protection: a comparison of the Boyd Orr Cohort of elderly subjects with a younger population sample. *Age Ageing* **36**, 521-526, doi:10.1093/ageing/afm107 (2007).
- 90 Langie, S. A. *et al.* Modulation of nucleotide excision repair in human lymphocytes by genetic and dietary factors. *Br J Nutr* **103**, 490-501, doi:10.1017/S0007114509992066 (2010).
- 91 Al-Serori, H. *et al.* Mobile phone specific electromagnetic fields induce transient DNA damage and nucleotide excision repair in serum-deprived human glioblastoma cells. *PLoS One* **13**, e0193677, doi:10.1371/journal.pone.0193677 (2018).
- 92 Soares, J. P. *et al.* Effects of combined physical exercise training on DNA damage and repair capacity: role of oxidative stress changes. *Age (Dordr)* **37**, 9799, doi:10.1007/s11357-015-9799-4 (2015).
- 93 Azqueta, A. *et al.* A comparative performance test of standard, medium- and high-throughput comet assays. *Toxicol In Vitro* **27**, 768-773, doi:10.1016/j.tiv.2012.12.006 (2013).
- 94 Collins, A. R. *et al.* The comet assay: topical issues. *Mutagenesis* **23**, 143-151, doi:10.1093/mutage/gem051 (2008).
- 95 Forchhammer, L. *et al.* Variation in the measurement of DNA damage by comet assay measured by the ECVAG inter-laboratory validation trial. *Mutagenesis* **25**, 113-123, doi:10.1093/mutage/gep048 (2010).
- 96 Azqueta, A., Langie, S. A., Slysokova, J. & Collins, A. R. Measurement of DNA base and nucleotide excision repair activities in mammalian cells and tissues using the comet assay--a methodological overview. *DNA Repair (Amst)* **12**, 1007-1010, doi:10.1016/j.dnarep.2013.07.011 (2013).

478 97 Gungor, N. *et al.* Lung inflammation is associated with reduced pulmonary nucleotide
479 excision repair in vivo. *Mutagenesis* **25**, 77-82, doi:10.1093/mutage/geb049 (2010).
480

Addition of
SAMPLE EXTRACT

**Recognition of
damaged base**

SUBSTRATE DNA

Nucleoids containing 8-oxoG lesions

INCUBATION REACTION

Incubation at 37°C

Base excision

AP site


DNA incision

**Denaturation during
alkaline treatment**

Electrophoresis

**COMET FORMATION
&
VISUALIZATION**


DNA migration into the comet tail

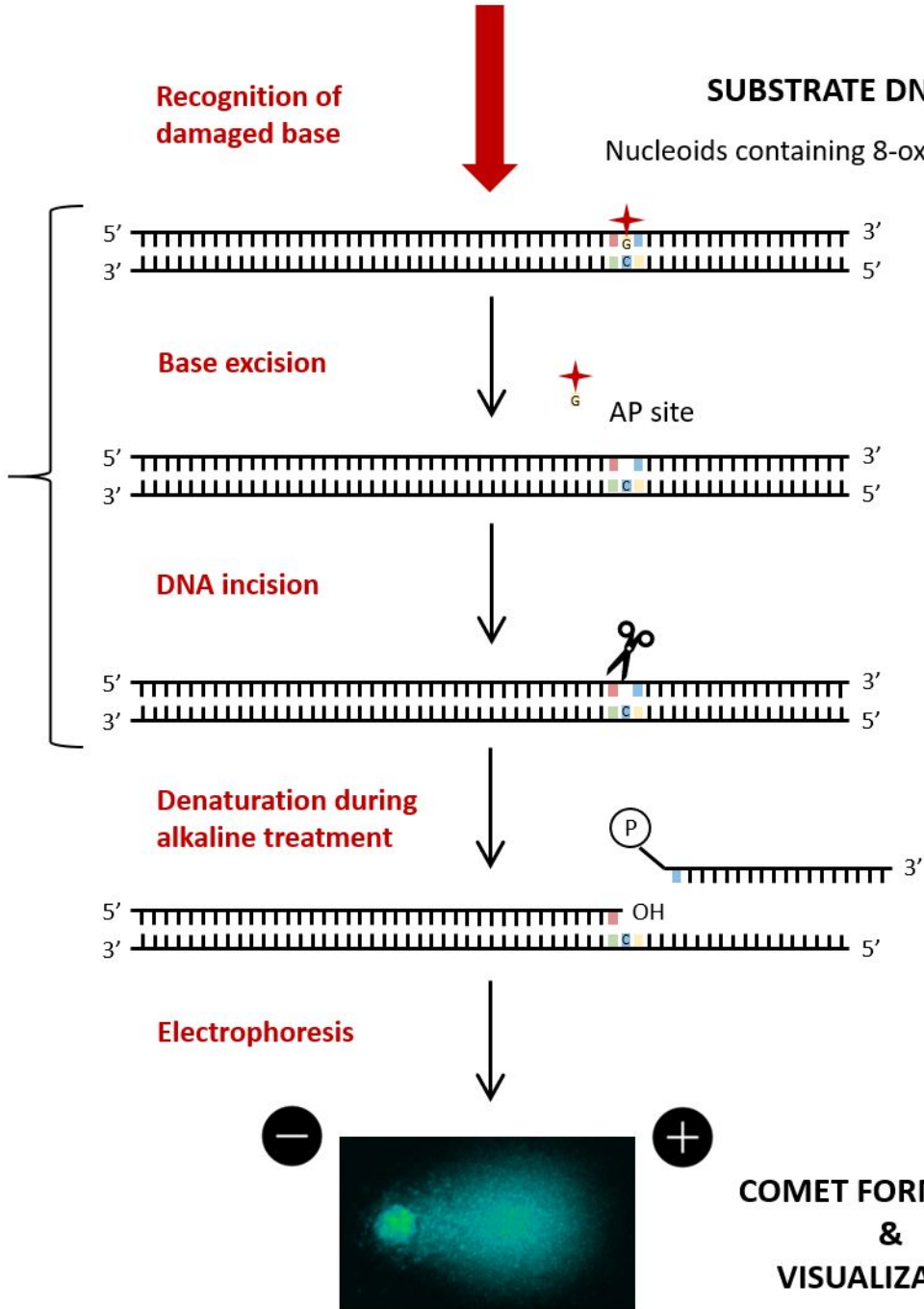


1

PREPARATION OF SUBSTRATE CELLS

Non-exposed cells

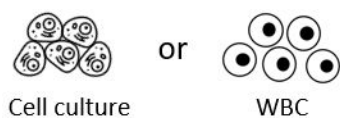


Exposed cells

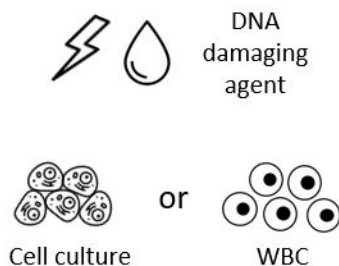

CHECKPOINT 1:
Assess levels of lesions in substrate cells

2

EXTRACT PREPARATION

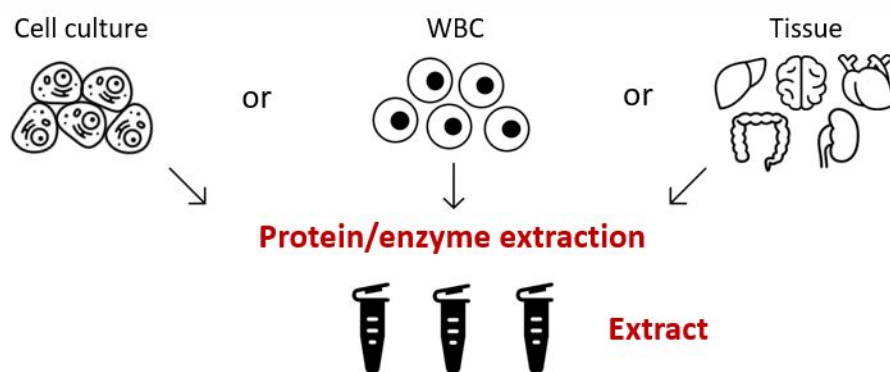


CHECKPOINT 2:
Measure protein concentration


3

EMBEDDING SUBSTRATE CELLS in LMP AGAROSE GEL & CELL LYSIS

Prepare 2 slides per sample and the buffer control

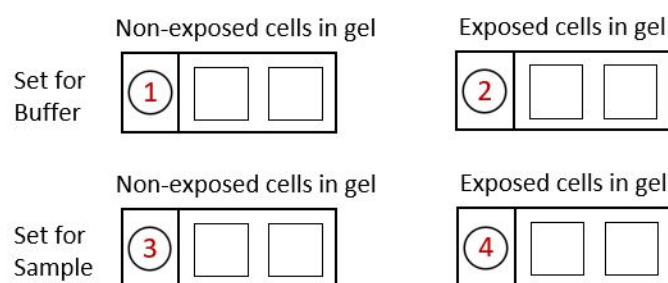


+ 2 extra slides
in each experiment
to incubate with
lesion-specific enzyme

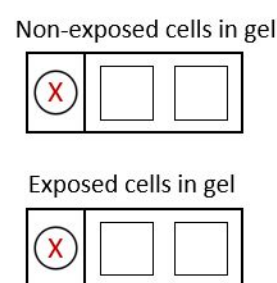

Immersion
of cells in
lysis
solution


Cell lysis

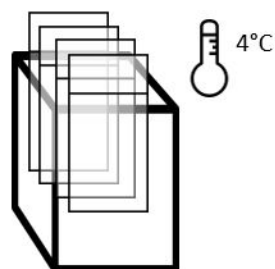


4

INCUBATION REACTION

Detection of DNA incision activity of the extract

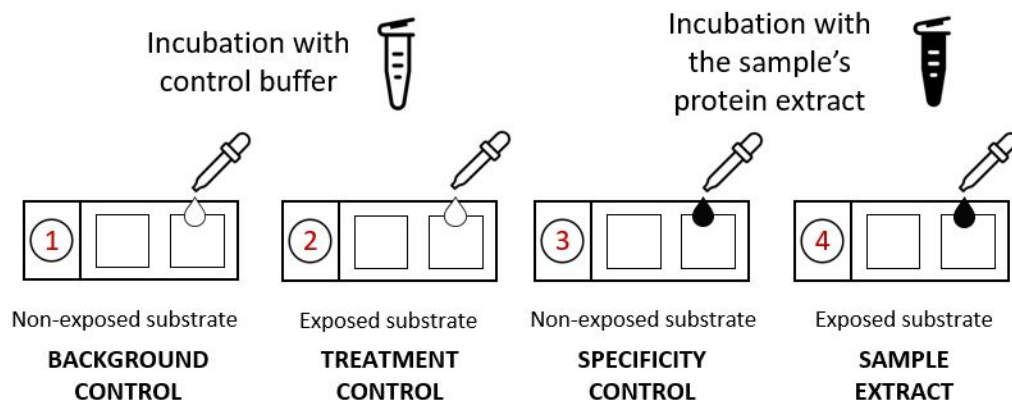

Incubation with
lesion-specific enzyme

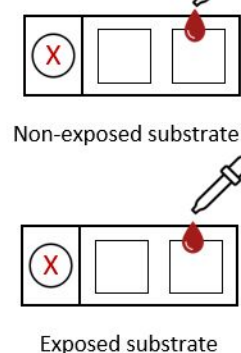

Incubation




Optional: Detection of DNA synthesis/ligation activity of the extract
Addition of dNTPs to the extract (parallel incubation)

INCUBATION REACTION
CONTROL

5

COMET FORMATION


Denaturation & Electrophoresis

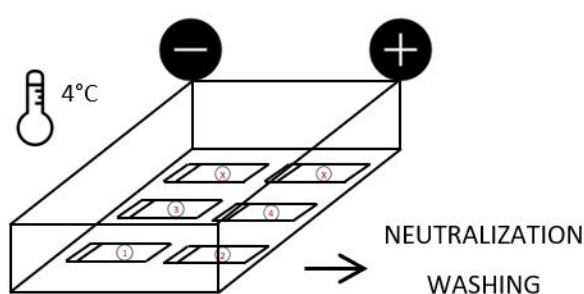

6

COMET VISUALIZATION & ANALYSIS

STAINING
&
MICROSCOPIC
EVALUATION


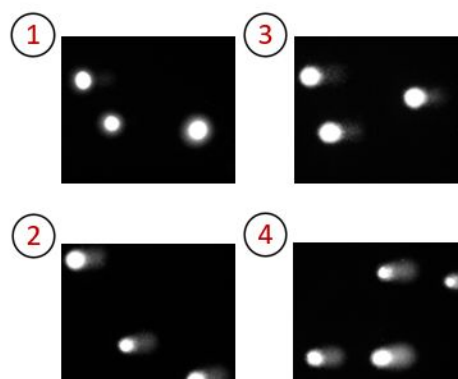


RESULTS

**Formula for calculating
repair-related DNA
incision activity:**


First subtract background value,
i.e. ① from all other mean values

DNA repair incision activity =
net TI ④ – net TI ③ – net TI ②

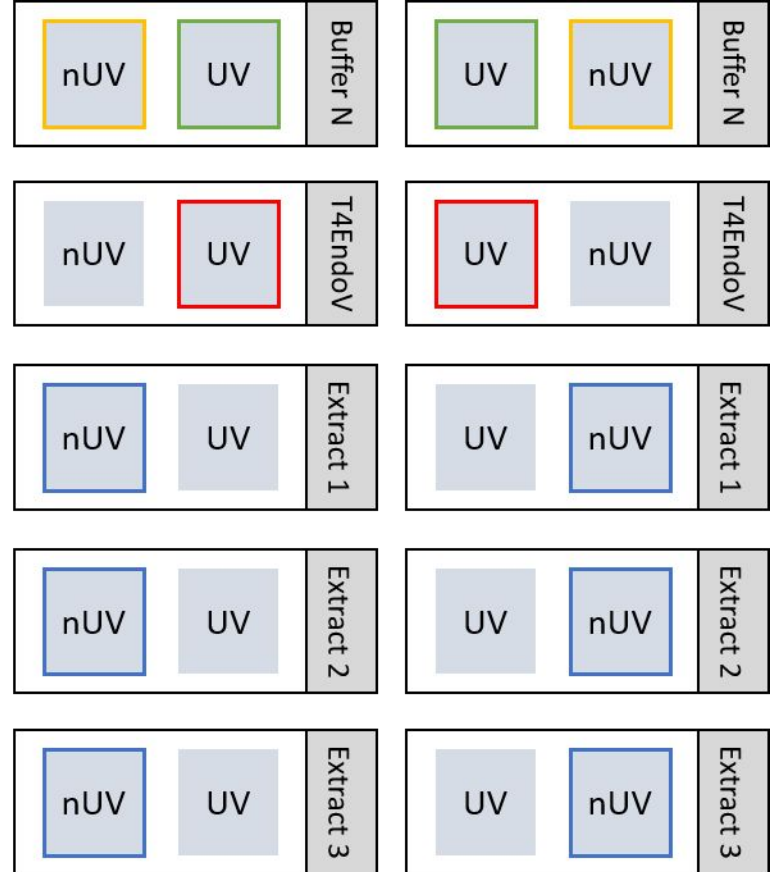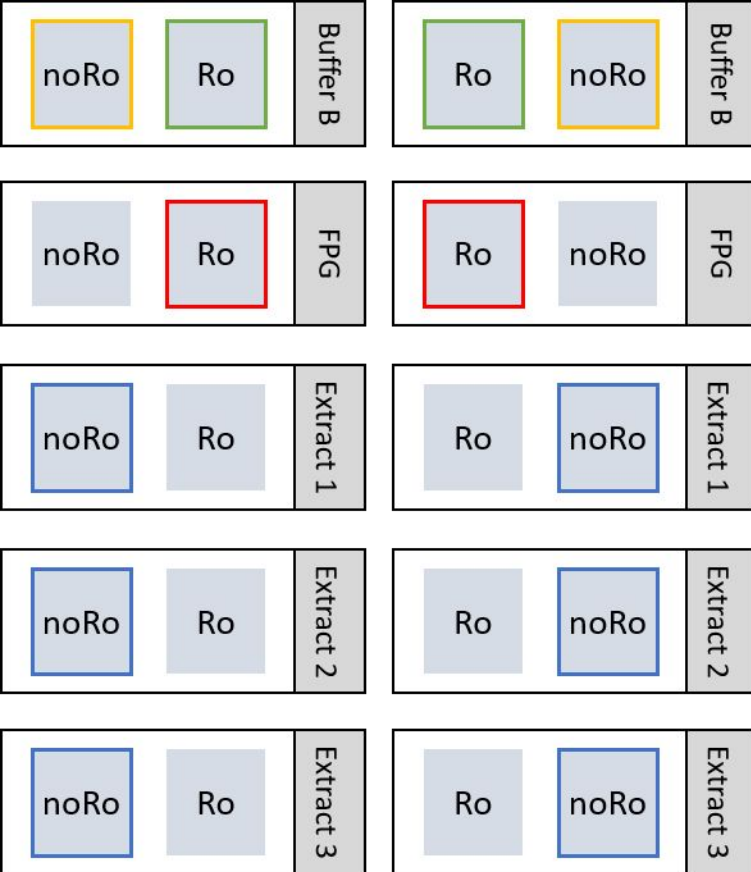


Legend:

- Background control (Yellow border)
- Treatment control (Green border)
- Specificity control (Blue border)
- Incubation reaction control (Red border)

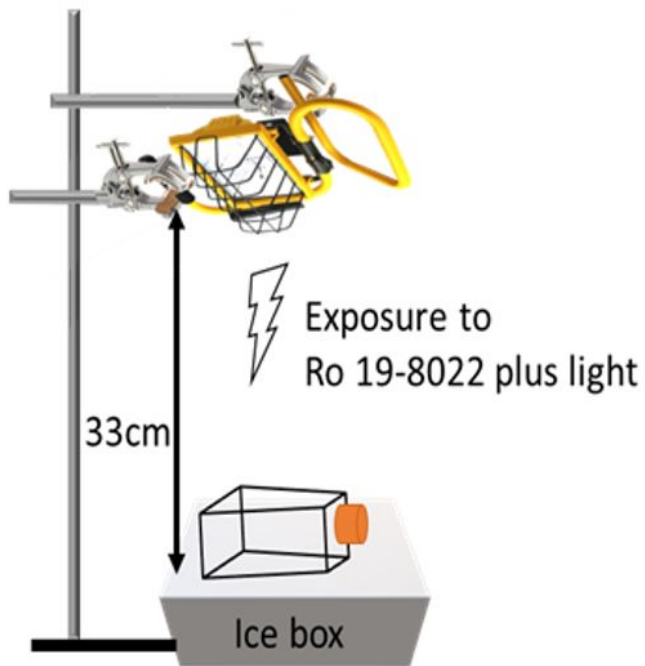



TOP

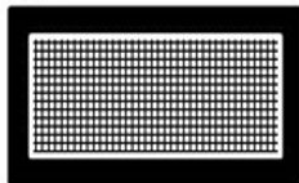

FRONT



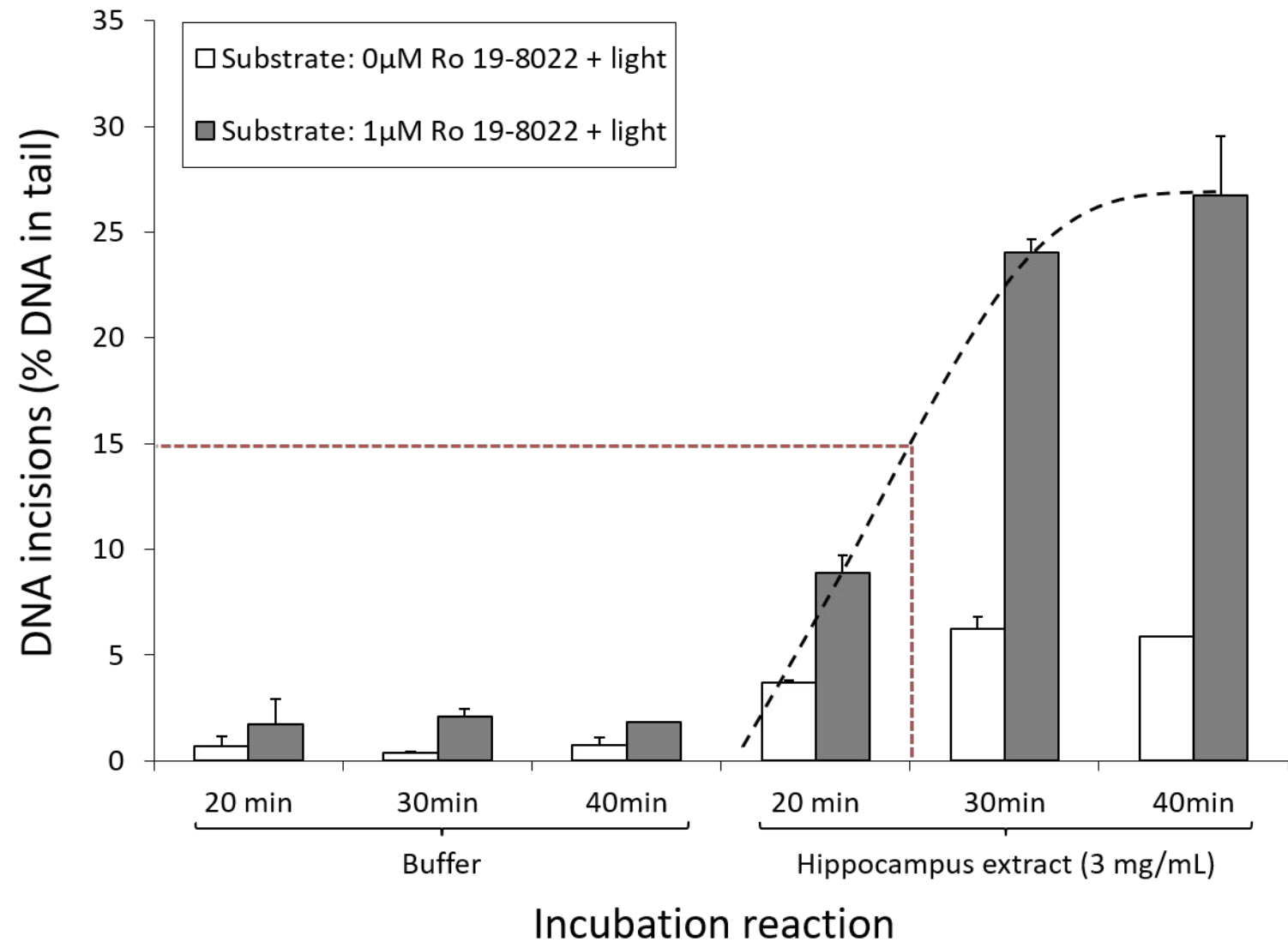


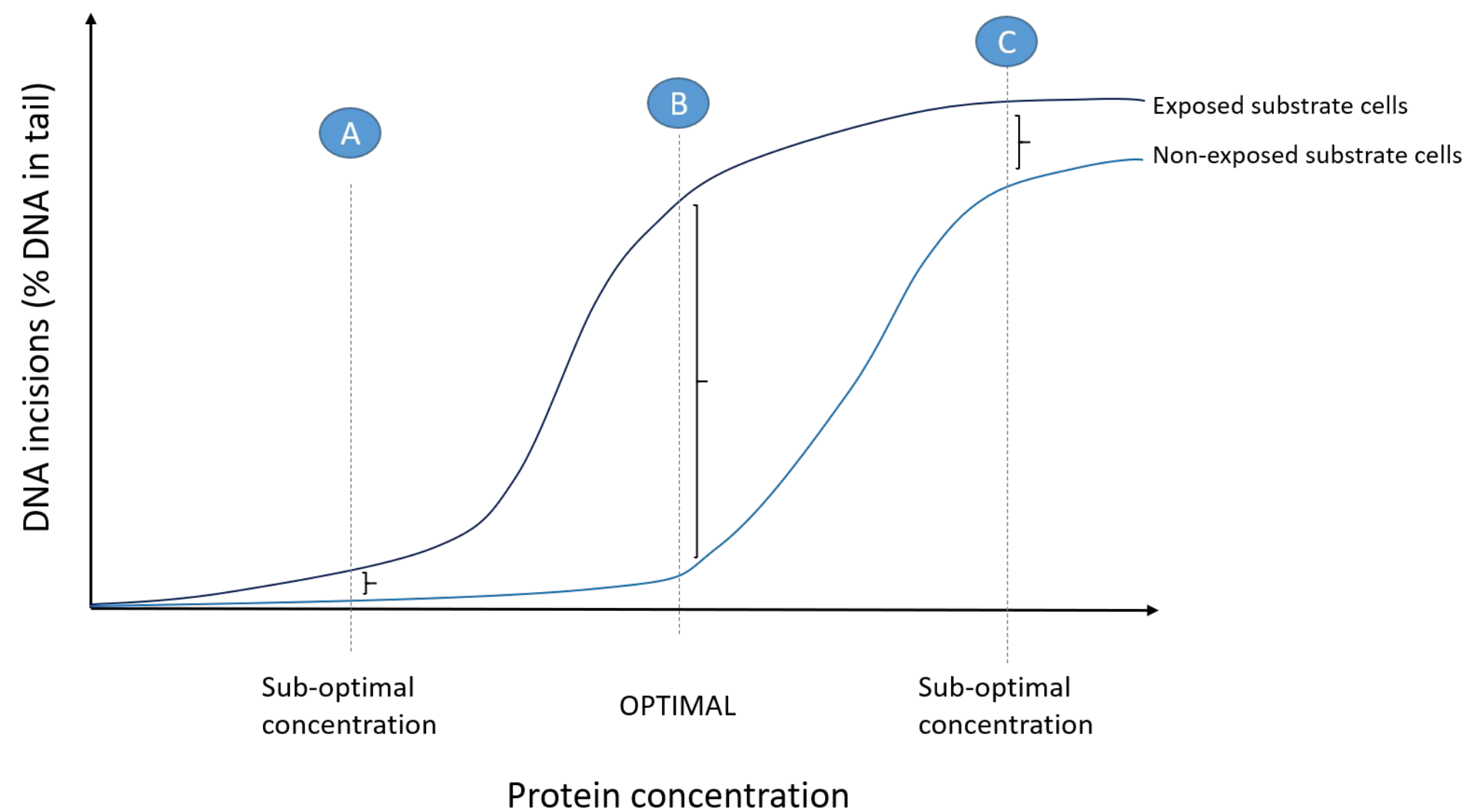


A)



B)


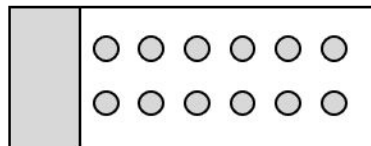


B 1 2 3 4 IC

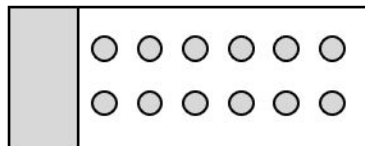

5 6 7 8 9 10

← Non-exposed substrate

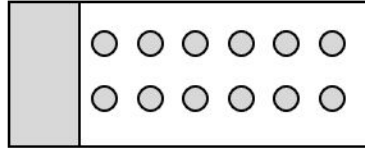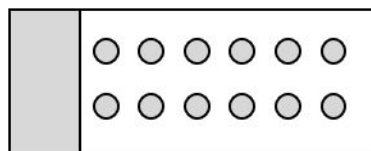


← Exposed substrate

Incubation with
control buffer




Incubation with
the sample's
protein extract

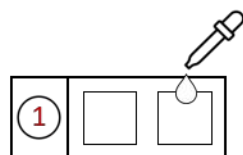


Non-exposed substrate

**BACKGROUND
CONTROL**

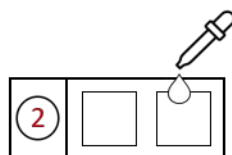


Exposed substrate

**TREATMENT
CONTROL**

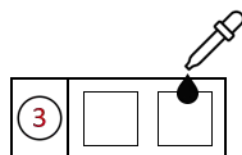


Non-exposed substrate



**SPECIFICITY
CONTROL**


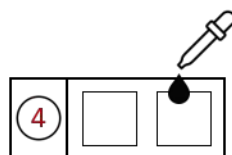


Exposed substrate


**SAMPLE
EXTRACT**

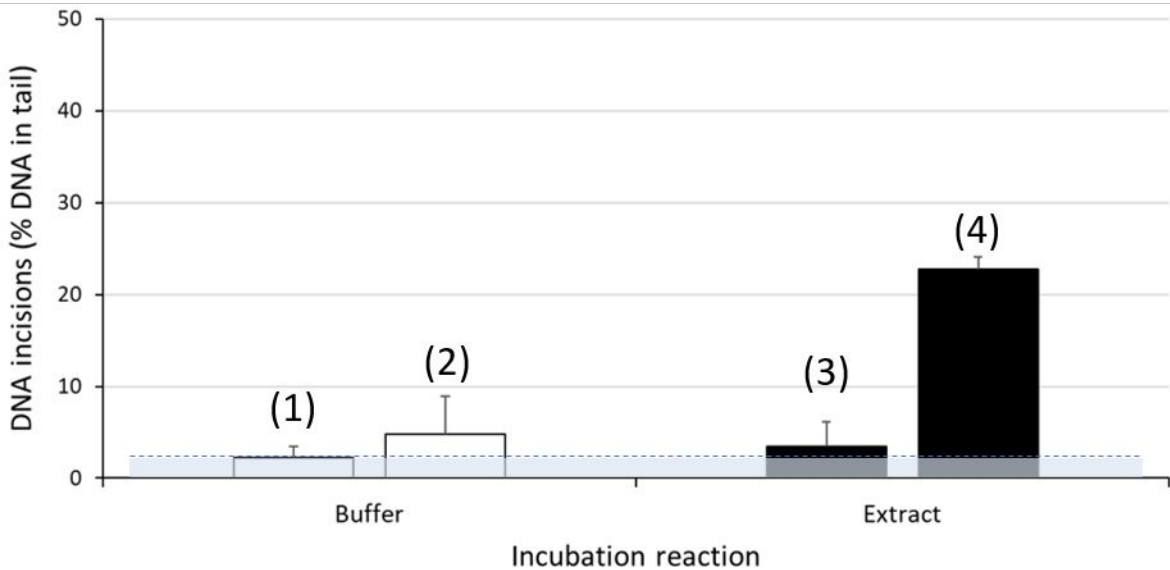

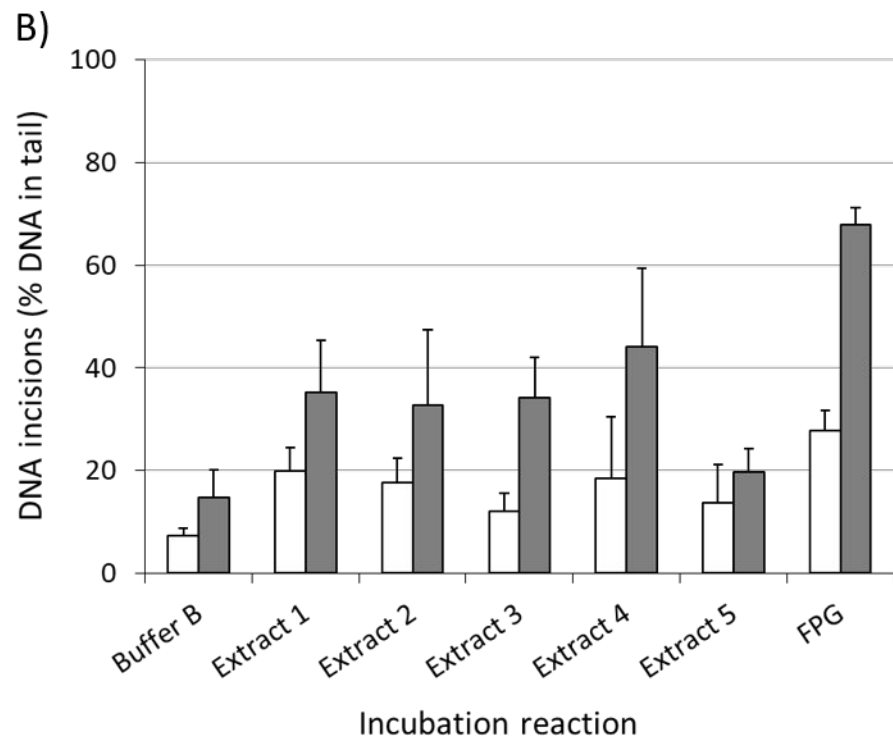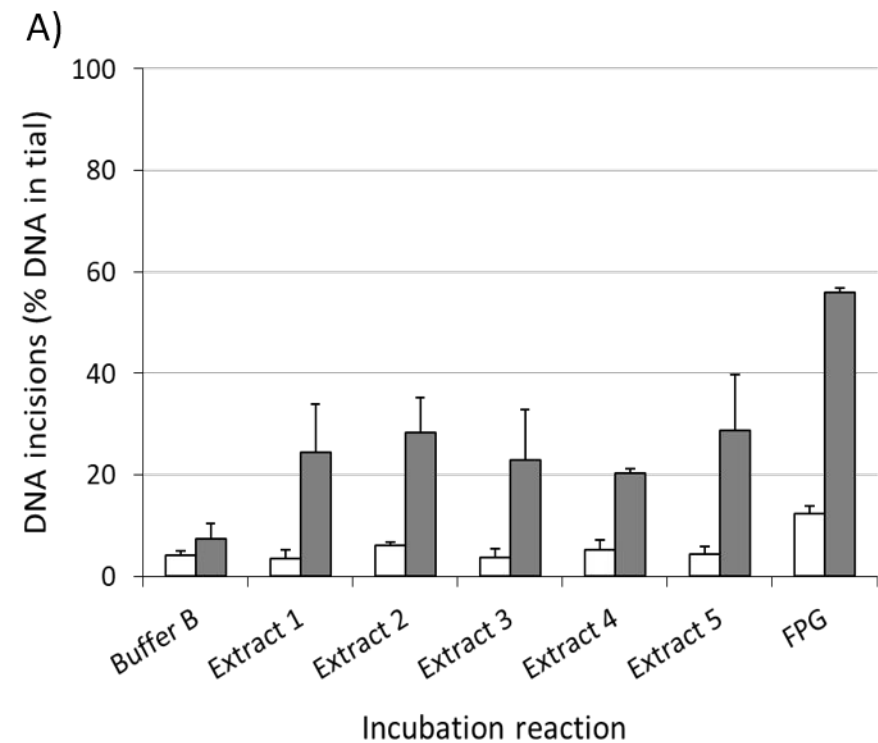

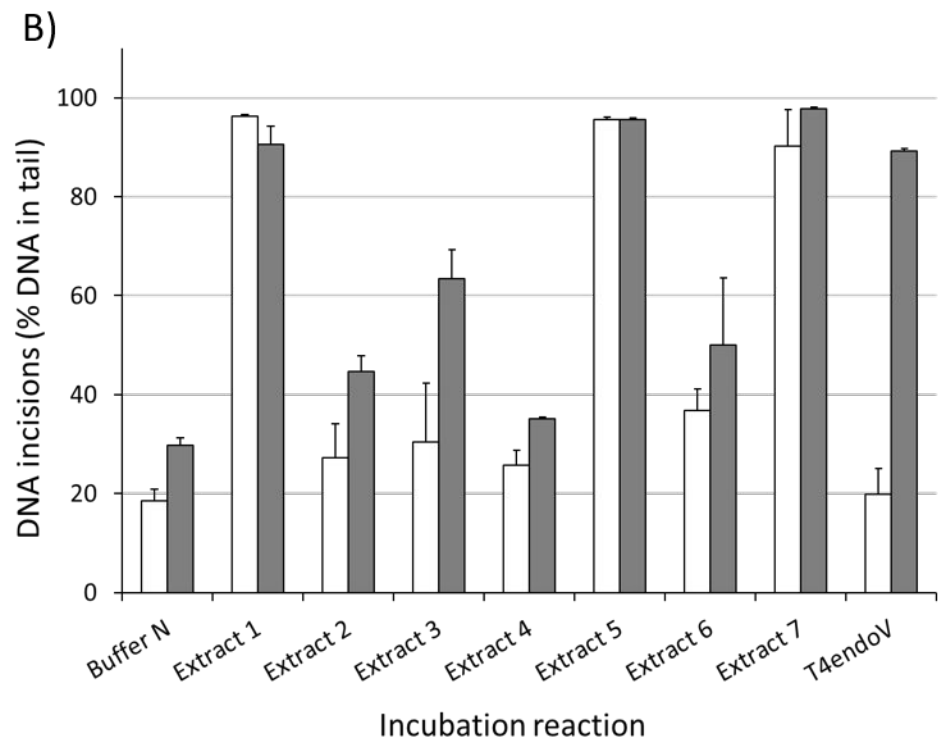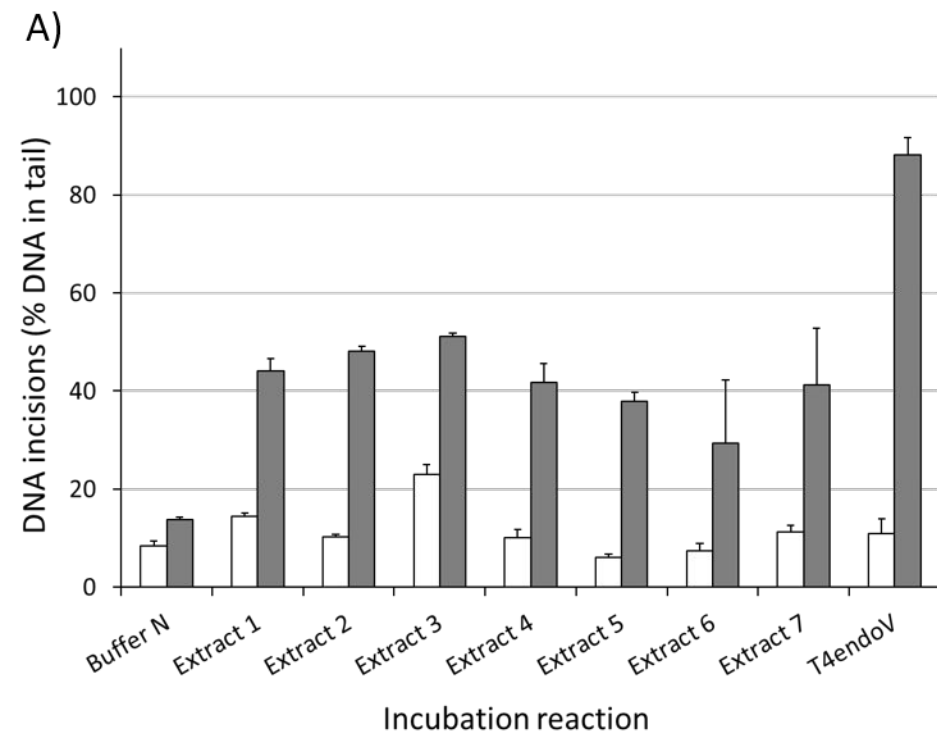

A)

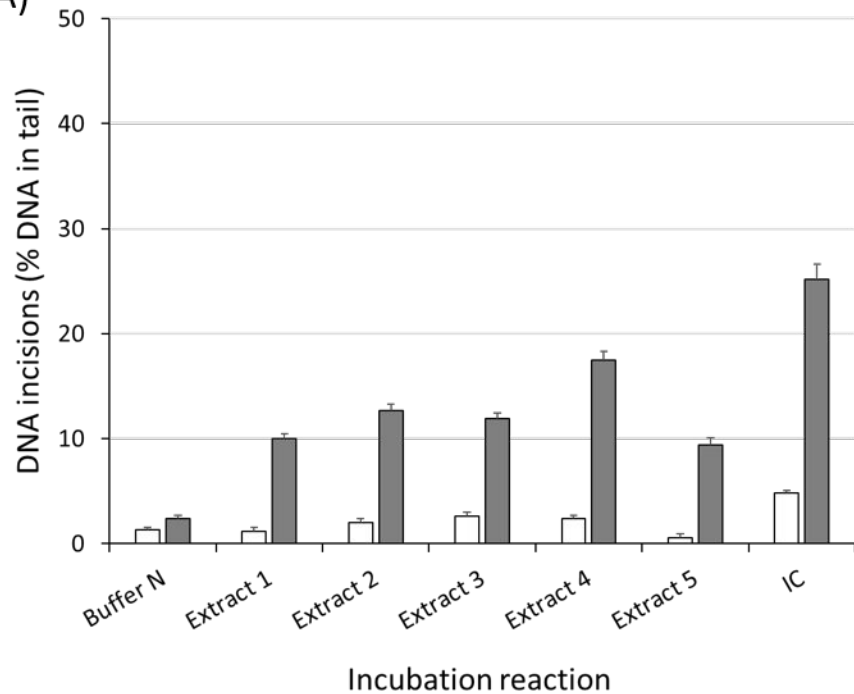

B)

